WDHD1 Rabbit pAb

Catalog No.: A15396



Basic Information

Observed MW

135kDa

Calculated MW

126kDa

Category

Primary antibody

Applications

ELISA,WB,IHC-P,IF/ICC

Cross-Reactivity

Human, Rat

Background

The protein encoded by this gene contains multiple N-terminal WD40 domains and a C-terminal high mobility group (HMG) box. WD40 domains are found in a variety of eukaryotic proteins and may function as adaptor/regulatory modules in signal transduction, pre-mRNA processing and cytoskeleton assembly. HMG boxes are found in many eukaryotic proteins involved in chromatin assembly, transcription and replication. Alternative splicing results in two transcript variants encoding different isoforms.

Recommended Dilutions

WB	1:200 - 1:2000
IHC-P	1:50 - 1:200
IF/ICC	1:50 - 1:200

Immunogen Information

Gene ID	Swiss Prot
11169	075717

Immunogen

Recombinant fusion protein containing a sequence corresponding to amino acids 830-1129 of human WDHD1 (NP_009017.1).

Synonyms

AND1; CTF4; AND-1; CHTF4; WDHD1

Contact

8		400-999-6126
\bowtie	Τ	cn.market@abclonal.com.cn
•	T	www.abclonal.com.cn


Product Information

SourceIsotypePurificationRabbitIgGAffinity purification

Storage

Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS with 0.01% thimerosal,50% glycerol,pH7.3.

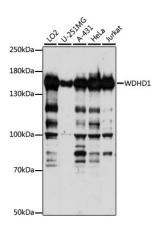


Western blot analysis of various lysates using WDHD1 Rabbit pAb (A15396) at 1:1000 dilution. Secondary antibody: HRP Goat Anti-Rabbit $\lg G$ (H+L) (AS014) at 1:10000 dilution.

Lysates/proteins: 25µg per lane.


Blocking buffer: 3% nonfat dry milk in TBST.

Detection: ECL Basic Kit (RM00020).

Exposure time: 10s.

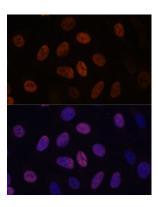

Immunofluorescence analysis of U-2 OS cells using WDHD1 Rabbit pAb (A15396) at dilution of 1:100. Secondary antibody: Cy3 Goat Anti-Rabbit IgG (H+L) (AS007) at 1:500 dilution. Blue: DAPI for nuclear staining.

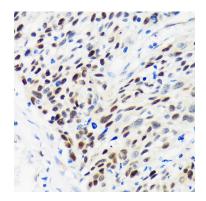

Immunohistochemistry analysis of paraffinembedded Human esophageal cancer using WDHD1 Rabbit pAb (A15396) at dilution of 1:100 (40x lens).Perform microwave antigen retrieval with 10 mM PBS buffer pH 7.2 before commencing with IHC staining protocol.