# **OGDHL Rabbit pAb**

Catalog No.: A15475 1 Publications



### **Basic Information**

#### **Observed MW**

120kDa

### **Calculated MW**

114kDa

### Category

Primary antibody

### **Applications**

ELISA,WB,IHC-P

#### **Cross-Reactivity**

Human, Mouse, Rat

# **Background**

The protein encoded by this gene is similar to oxoglutarate dehydrogenase (OGDH) of the OGDH complex, which degrades glucose and glutamate. This gene encodes several isoforms, including some that appear to localize to mitochondria. The encoded protein down-regulates the AKT signaling cascade and can suppress the growth of cervical cancer cells.

# **Recommended Dilutions**

**WB** 1:200 - 1:2000

**IHC-P** 1:50 - 1:200

# **Immunogen Information**

**Gene ID**Swiss Prot
55753
Q9ULD0

#### **Immunogen**

Recombinant fusion protein containing a sequence corresponding to amino acids 1-110 of human OGDHL ( $NP_060715.2$ ).

# **Synonyms**

YOBELN; OGDHL

### **Contact**

<b>a</b>		400-999-6126
$\bowtie$		cn.market@abclonal.com.cn
$\Box$	ī	www.abclonal.com.cn


### **Product Information**

SourceIsotypePurificationRabbitIgGAffinity purification

### **Storage**

Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS with 0.01% thimerosal,50% glycerol,pH7.3.

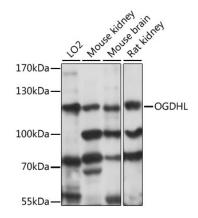


Western blot analysis of various lysates using OGDHL Rabbit pAb (A15475) at 1:1000 dilution. Secondary antibody: HRP Goat Anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution. Lysates/proteins:  $25\mu g$  per lane.


Blocking buffer: 3% nonfat dry milk in TBST.

Detection: ECL Basic Kit (RM00020).

Exposure time: 5s.

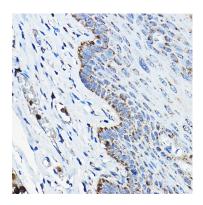

Immunohistochemistry analysis of paraffinembedded human esophageal using OGDHL Rabbit pAb (A15475) at dilution of 1:200 (40x lens).Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.

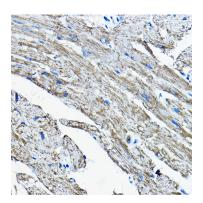

Immunohistochemistry analysis of paraffinembedded mouse heart using OGDHL Rabbit pAb (A15475) at dilution of 1:200 (40x lens).Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.