Leader in Biomolecular Solutions for Life Science



BRD2 Rabbit pAb

Catalog No.: A16241 2 Publications

Basic Information

Observed MW 110kDa

Calculated MW 88kDa

Category Primary antibody

Applications ELISA,WB,IHC-P

Cross-Reactivity Human, Mouse, Rat

Background

This gene encodes a transcriptional regulator that belongs to the BET (bromodomains and extra terminal domain) family of proteins. This protein associates with transcription complexes and with acetylated chromatin during mitosis, and it selectively binds to the acetylated lysine-12 residue of histone H4 via its two bromodomains. The gene maps to the major histocompatability complex (MHC) class II region on chromosome 6p21.3, but sequence comparison suggests that the protein is not involved in the immune response. This gene has been implicated in juvenile myoclonic epilepsy, a common form of epilepsy that becomes apparent in adolescence. Multiple alternatively spliced variants have been described for this gene.

Recommended Dilutions

Immunogen Information

WB	1:500 - 1:1000	Gene ID	Swiss Prot	
IHC-P	1:50 - 1:200	6046	P25440	
		I		

Immunogen

A synthetic peptide corresponding to a sequence within amino acids 301-400 of human BRD2 (NP_005095.1).

Synonyms

FSH; NAT; RNF3; FSRG1; RING3; FSHRG1; D6S113E; O27.1.1; BRD2-IT1; BRD2


a 400-999-6126 x cn.market@abclonal.com.cn y www.abclonal.com.cn


Product Information

Source Rabbit **Isotype** IgG **Purification** Affinity purification


Storage

Store at -20°C. Avoid freeze / thaw cycles. Buffer: PBS with 0.05% proclin300,50% glycerol,pH7.3.

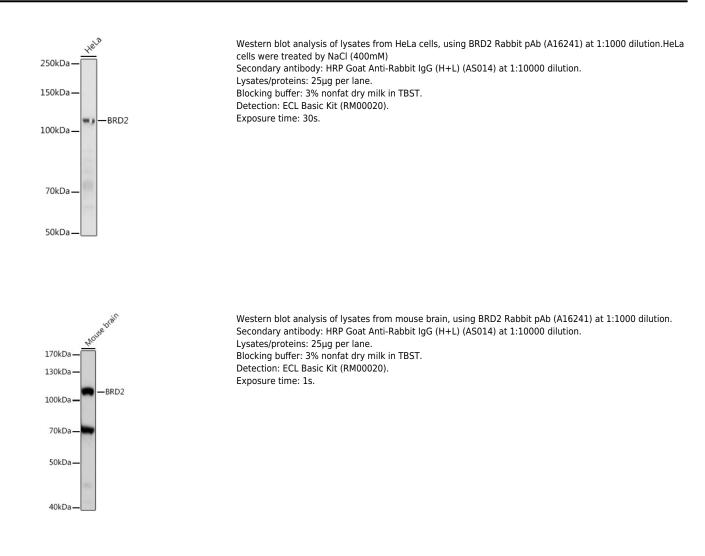

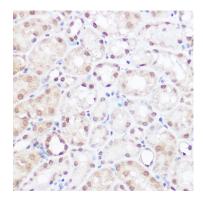

Immunohistochemistry analysis of paraffinembedded mouse kidney using BRD2 Rabbit pAb (A16241) at dilution of 1:50 (40x lens).Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.

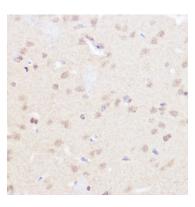

Immunohistochemistry analysis of paraffinembedded rat brain using BRD2 Rabbit pAb (A16241) at dilution of 1:50 (40x lens).Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.