METTL3 Rabbit pAb

Catalog No.: A21572 1 Publications



Basic Information

Observed MW

75kDa

Calculated MW

64kDa

Category

Primary antibody

Applications

ELISA,WB,IHC-P,IP

Cross-Reactivity

Human, Mouse, Rat

Background

This gene encodes the 70 kDa subunit of MT-A which is part of N6-adenosine-methyltransferase. This enzyme is involved in the posttranscriptional methylation of internal adenosine residues in eukaryotic mRNAs, forming N6-methyladenosine.

Recommended Dilutions

WB	1:500 - 1:1000
IHC-P	1:10 - 1:20
IHC-P	1:500 - 1:1000

Immunogen Information

Gene ID	Swiss Prot
56339	Q86U44

Immunogen

Recombinant fusion protein containing a sequence corresponding to amino acids 1-156 of human METTL3 (NP_062826.2).

Synonyms

M6A; IME4; Spo8; MT-A70; hMETTL3; METTL3

Contact

6	400-999-6126
\bowtie	cn.market@abclonal.com.cn
•	www.abclonal.com.cn


Product Information

SourceIsotypePurificationRabbitIgGAffinity purification

Storage

Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS with 0.05% proclin300,50% glycerol,pH7.3.

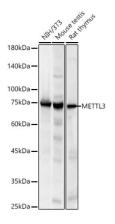


Western blot analysis of various lysates, using METTL3 Rabbit pAb (A21572) at 1:1000 dilution. Secondary antibody: HRP Goat Anti-Rabbit lgG (H+L) (AS014) at 1:10000 dilution.

Lysates/proteins: 25µg per lane.


Blocking buffer: 3% nonfat dry milk in TBST.

Detection: ECL Basic Kit (RM00020).


Exposure time: 20s.

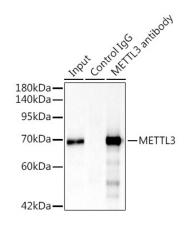

Immunoprecipitation analysis of 300 μ g extracts of 293T cells using 3 μ g METTL3 Rabbit pAb (A21572). Western blot was performed from the immunoprecipitate using METTL3 Rabbit pAb (A21572) at a dilution of 1,1000.

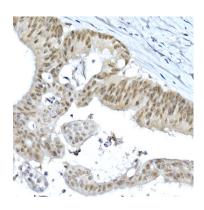

Immunohistochemistry analysis of paraffinembedded Human colon carcinoma using METTL3 Rabbit pAb (A21572) at dilution of 1:20 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate Bufferr (pH 6.0) prior to IHC staining.

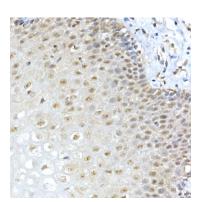

Immunohistochemistry analysis of paraffinembedded Human esophagus using METTL3 Rabbit pAb (A21572) at dilution of 1:20 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate Bufferr (pH 6.0) prior to IHC staining.