NEDD4 Rabbit pAb

Catalog No.: A0552 3 Publications



Basic Information

Observed MW

120kDa

Calculated MW

149kDa

Category

Primary antibody

Applications

ELISA,WB,IHC-P

Cross-Reactivity

Human

Background

This gene is the founding member of the NEDD4 family of HECT ubiquitin ligases that function in the ubiquitin proteasome system of protein degradation. The encoded protein contains an N-terminal calcium and phospholipid binding C2 domain followed by multiple tryptophan-rich WW domains and, a C-terminal HECT ubiquitin ligase catalytic domain. It plays critical role in the regulation of a number of membrane receptors, endocytic machinery components and the tumor suppressor PTEN.

Recommended Dilutions

WB 1:500 - 1:2000

IHC-P 1:50 - 1:200

Immunogen Information

Gene ID4734

Swiss Prot
P46934

Immunogen

Recombinant fusion protein containing a sequence corresponding to amino acids 160-460 of human NEDD4 (NP_940682.2).

Synonyms

RPF1; NEDD4-1; NEDD4

Contact

a	400-999-6126
×	cn.market@abclonal.com.cn
\overline{a}	www.ahclonal.com.cn


Product Information

SourceIsotypePurificationRabbitIgGAffinity purification

Storage

Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS with 0.02% sodium azide,50% glycerol,pH7.3.

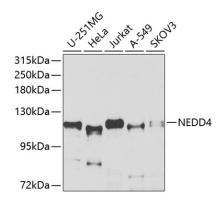


Western blot analysis of extracts of various cell lines, using NEDD4 antibody (A0552) at 1:1000 dilution. Secondary antibody: HRP Goat Anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.


Lysates/proteins: 25µg per lane.

Blocking buffer: 3% nonfat dry milk in TBST. Detection: ECL Enhanced Kit (RM00021).

Exposure time: 30s.

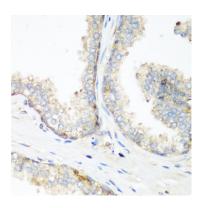

Immunohistochemistry analysis of paraffinembedded human prostate using NEDD4 antibody (A0552) at dilution of 1:200 (40x lens).Perform microwave antigen retrieval with 10 mM PBS buffer pH 7.2 before commencing with IHC staining protocol.

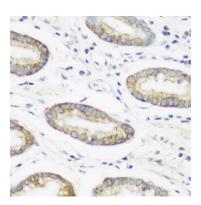

Immunohistochemistry analysis of paraffinembedded human gastric using NEDD4 antibody (A0552) at dilution of 1:200 (40x lens).Perform microwave antigen retrieval with 10 mM PBS buffer pH 7.2 before commencing with IHC staining protocol.