# LRRK2 Rabbit pAb

Catalog No.: A0859



## **Basic Information**

#### **Observed MW**

290kDa

### **Calculated MW**

286kDa

### Category

Primary antibody

### **Applications**

ELISA,WB,IHC-P

### **Cross-Reactivity**

Human, Mouse, Rat

## **Background**

This gene is a member of the leucine-rich repeat kinase family and encodes a protein with an ankryin repeat region, a leucine-rich repeat (LRR) domain, a kinase domain, a DFG-like motif, a RAS domain, a GTPase domain, a MLK-like domain, and a WD40 domain. The protein is present largely in the cytoplasm but also associates with the mitochondrial outer membrane. Mutations in this gene have been associated with Parkinson disease-8.

## **Recommended Dilutions**

**WB** 1:500 - 1:1000

**IHC-P** 1:50 - 1:200

## **Immunogen Information**

 Gene ID
 Swiss Prot

 120892
 Q5S007

### **Immunogen**

Recombinant fusion protein containing a sequence corresponding to amino acids 750-910 of human LRRK2 (NP\_940980.3).

### **Synonyms**

PARK8; RIPK7; ROCO2; AURA17; DARDARIN; LRRK2

## **Contact**

<u>a</u>		400-999-6126
$\bowtie$		cn.market@abclonal.com.cn
<u> </u>	Π	www.ahclonal.com.cn


## **Product Information**

SourceIsotypePurificationRabbitIgGAffinity purification

## Storage

Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS with 0.02% sodium azide,50% glycerol,pH7.3.

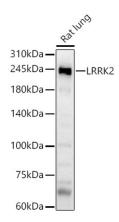


Western blot analysis of Rat lung, using LRRK2 antibody (A0859) at 1:900 dilution. Secondary antibody: HRP Goat Anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution. Lysates/proteins:  $25\mu g$  per lane.


Blocking buffer: 3% nonfat dry milk in TBST.

Detection: ECL Basic Kit (RM00020).

Exposure time: 60s.

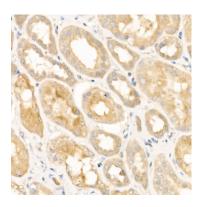

Immunohistochemistry analysis of paraffinembedded human kidney using LRRK2 Rabbit pAb (A0859) at dilution of 1:20 (40x lens).Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.

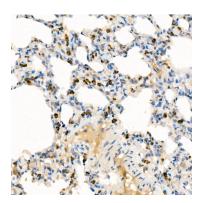

Immunohistochemistry analysis of paraffinembedded rat lung using LRRK2 Rabbit pAb (A0859) at dilution of 1:20 (40x lens).Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.