# **ILK Rabbit pAb**

Catalog No.: A0901



## **Basic Information**

#### **Observed MW**

48kDa

### **Calculated MW**

51kDa

### Category

Primary antibody

### **Applications**

ELISA,WB,IHC-P

### **Cross-Reactivity**

Human, Mouse, Rat

## **Background**

This gene encodes a protein with a kinase-like domain and four ankyrin-like repeats. The encoded protein associates at the cell membrane with the cytoplasmic domain of beta integrins, where it regulates integrin-mediated signal transduction. Activity of this protein is important in the epithelial to mesenchymal transition, and over-expression of this gene is implicated in tumor growth and metastasis. Alternative splicing results in multiple transcript variants.

## **Recommended Dilutions**

**WB** 1:500 - 1:1000

**IHC-P** 1:50 - 1:100

## **Immunogen Information**

**Gene ID**Swiss Prot
3611
Q13418

### **Immunogen**

Recombinant fusion protein containing a sequence corresponding to amino acids 1-260 of human ILK (NP\_001014795.1).

### **Synonyms**

P59; ILK-1; ILK-2; p59ILK; HEL-S-28; ILK

## **Contact**

2		400-999-6126
$\bowtie$		cn.market@abclonal.com.cn
•	T	www.abclonal.com.cn


## **Product Information**

SourceIsotypePurificationRabbitIgGAffinity purification

### **Storage**


Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS with 0.02% sodium azide,50% glycerol,pH7.3.

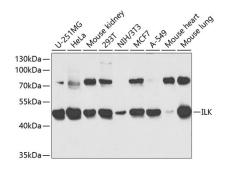

Western blot analysis of various lysates using ILK Rabbit pAb (A0901) at 1:1000 dilution. Secondary antibody: HRP Goat Anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution. Lysates/proteins: 25µg per lane.

Blocking buffer: 3% nonfat dry milk in TBST.

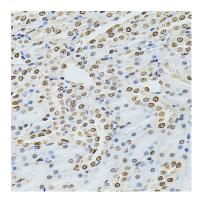

Immunohistochemistry analysis of paraffinembedded rat kidney using ILK Rabbit pAb (A0901) at dilution of 1:200 (40x lens).Perform microwave antigen retrieval with 10 mM PBS buffer pH 7.2 before commencing with IHC staining protocol.



Immunohistochemistry analysis of paraffinembedded mouse brain using ILK Rabbit pAb (A0901) at dilution of 1:200 (40x lens).Perform microwave antigen retrieval with 10 mM PBS buffer pH 7.2 before commencing with IHC staining protocol.