MTCO2 Rabbit pAb

Catalog No.: A11154 1 Publications



Basic Information

Observed MW

24kDa

Calculated MW

26kDa

Category

Primary antibody

Applications

WB,IHC-P,IF/ICC,(ELISA)

Cross-Reactivity

Human, Mouse, Rat

Background

Predicted to enable copper ion binding activity and oxidoreductase activity. Predicted to contribute to cytochrome-c oxidase activity. Predicted to be involved in ATP synthesis coupled electron transport; positive regulation of cellular biosynthetic process; and positive regulation of necrotic cell death. Located in mitochondrion. Is expressed in embryo; epiblast; heart; liver; and metanephros. Orthologous to several human genes including MTCO2P12 (MT-CO2 pseudogene 12).

Recommended Dilutions

WB 1:500 - 1:2000

IHC-P 1:100 - 1:500

IF/ICC 1:50 - 1:200

ELISA Recommended starting

concentration is 1 µg/mL.

Please optimize the
concentration based on
your specific assay
requirements.

Immunogen Information

Gene IDSwiss Prot
17709
P00403

Immunogen

A synthetic peptide corresponding to a sequence within amino acids 101-200 of mouse MTCO2 (NP_904331.1).

Synonyms

COX2; MTCO2

Contact

<u>a</u>		400-999-6126
\bowtie		cn.market@abclonal.com.cn
•	Ī	www.abclonal.com.cn


Product Information

SourceIsotypePurificationRabbitIgGAffinity purification


Storage

Store at -20°C. Avoid freeze / thaw cycles.


Buffer: PBS with 0.02% sodium azide,50% glycerol,pH7.3.

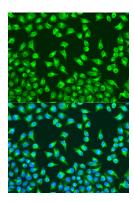

Immunofluorescence analysis of U2OS cells using MTCO2 Rabbit pAb (A11154) at dilution of 1:100. Secondary antibody: Cy3-conjugated Goat anti-Rabbit IgG (H+L) (AS007) at 1:500 dilution. Blue: DAPI for nuclear staining.

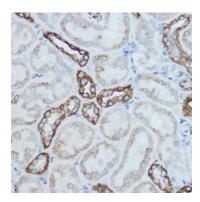

Immunohistochemistry analysis of paraffinembedded Mouse kidney using MTCO2 Rabbit pAb (A11154) at dilution of 1:500 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate Bufferr (pH 6.0) prior to IHC staining.

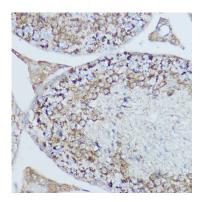

Immunohistochemistry analysis of paraffinembedded Rat testis using MTCO2 Rabbit pAb (A11154) at dilution of 1:500 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate Bufferr (pH 6.0) prior to IHC staining.