

p75NTR/NGFR Rabbit pAb

Catalog No.: A11169SP

Basic Information

Observed MW

60-75 kDa

Calculated MW

45 kDa/36 kDa

Category

Primary antibody

Applications

WB, IHC-P, ELISA

Cross-Reactivity

Human, Mouse, Rat

Background

Nerve growth factor receptor contains an extracellular domain containing four 40-amino acid repeats with 6 cysteine residues at conserved positions followed by a serine/threonine-rich region, a single transmembrane domain, and a 155-amino acid cytoplasmic domain. The cysteine-rich region contains the nerve growth factor binding domain.

Recommended Dilutions

WB 1:1000 - 1:2000

IHC-P 1:500 - 1:2000

ELISA Recommended starting concentration is 1 µg/mL. Please optimize the concentration based on your specific assay requirements. For high-ratio antibody dilutions ($\geq 1:10000$) a sequential dilution method is strongly recommended to ensure measurement accuracy.

Immunogen Information

Gene ID

4804/18053

Swiss Prot

P08138/Q9Z0W1

Immunogen

This information is considered to be commercially sensitive.

Synonyms

CD271; p75NTR; TNFRSF16; p75(NTR); Gp80-LNGFR

Product Information

Source

Rabbit

Isotype

IgG

Purification

Affinity purification

Storage

Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS, pH 7.3, containing 50% glycerol. Preserved with Proclin300 or sodium azide.

May contain 0.05% BSA as specified on the Certificate of Analysis.


Contact

☎ | 400-999-6126


✉ | cn.market@abclonal.com.cn

🌐 | www.abclonal.com.cn

Validation Data

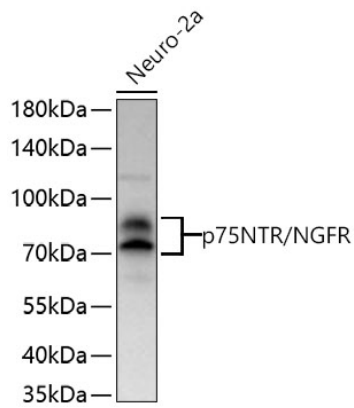

Western blot analysis of lysates from Neuro-2a cells using p75NTR/NGFR Rabbit pAb (A11169SP) at 1:1000 dilution incubated at room temperature for 1.5 hours.
Secondary antibody: HRP-conjugated Goat anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.
Lysates/proteins: 25 µg per lane.
Blocking buffer: 3% nonfat dry milk in TBST.
Detection: ECL Basic Kit (RM00020).
Exposure time: 45 s.

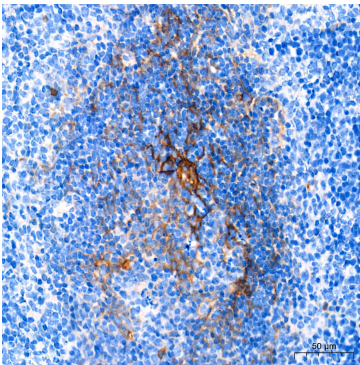

Immunohistochemistry analysis of paraffin-embedded Mouse spleen tissue using p75NTR/NGFR Rabbit pAb (A11169SP) at a dilution of 1:1500 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer (pH 9.0) prior to IHC staining.