[KD Validated] FAK Rabbit pAb

Catalog No.: A11195 21 Publications



Basic Information

Observed MW 125Kda/119kDa

Calculated MW 119kDa

Category Primary antibody

Applications WB,IHC-P,IF/ICC,IP,ELISA

Cross-Reactivity Human, Mouse, Rat

Background

This gene encodes a cytoplasmic protein tyrosine kinase which is found concentrated in the focal adhesions that form between cells growing in the presence of extracellular matrix constituents. The encoded protein is a member of the FAK subfamily of protein tyrosine kinases but lacks significant sequence similarity to kinases from other subfamilies. Activation of this gene may be an important early step in cell growth and intracellular signal transduction pathways triggered in response to certain neural peptides or to cell interactions with the extracellular matrix. Several transcript variants encoding different isoforms have been found for this gene.

Recommended Dilutions

| WB | 1:500 - 1:1000 |
|--------|--|
| IHC-P | 1:50 - 1:200 |
| IF/ICC | 1:50 - 1:200 |
| IP | 0.5µg-4µg antibody for 200µg-400µg extracts of whole cells |
| ELISA | Recommended starting concentration is 1 µg/mL. Please optimize the concentration based on your specific assay requirements. |

Contact

| 6 | 400-999-6126 |
|-----------|---------------------------|
| \bowtie | cn.market@abclonal.com.cn |
| € | www.abclonal.com.cn |

Immunogen Information

Gene ID 5747

Swiss Prot Q05397

Immunogen

Recombinant protein (or fragment). This information is considered to be commercially sensitive.

Synonyms

FAK; FADK; FAK1; FRNK; FADK 1; PPP1R71; p125FAK; pp125FAK

Product Information

Source Rabbit **Isotype** IgG Purification Affinity purification

Storage


Store at -20°C. Avoid freeze / thaw cycles. Buffer: PBS with 0.09% Sodium azide,50% glycerol,pH7.3.

150kDa

100kDa


70kDa

50kDa

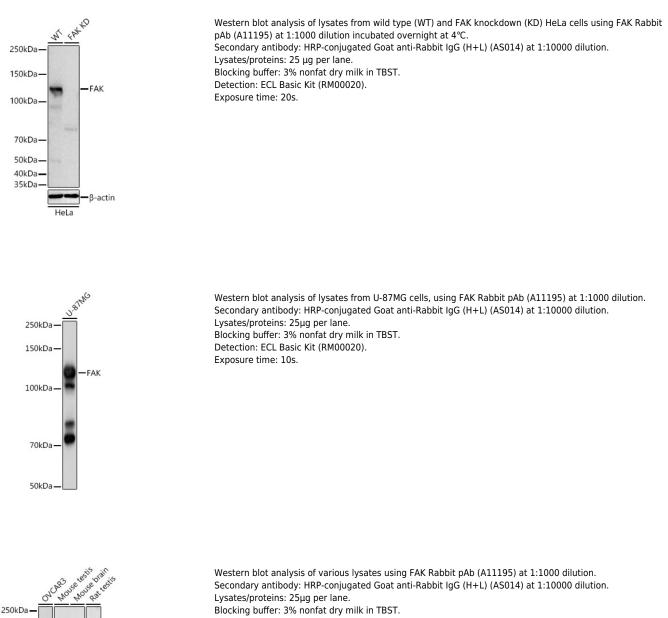


Blocking buffer: 3% nonfat dry milk i Detection: ECL Basic Kit (RM00020). Exposure time: 30s.


FAK

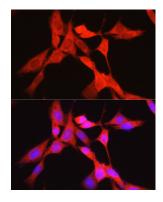


Immunofluorescence analysis of NIH/3T3 cells using FAK Rabbit pAb (A11195) at dilution of 1:50 (40x lens). Secondary antibody: Cy3-conjugated Goat anti-Rabbit IgG (H+L) (AS007) at 1:500 dilution. Blue: DAPI for nuclear staining.

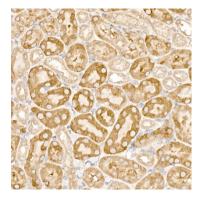


Immunohistochemistry analysis of paraffinembedded Mouse kidney using FAK Rabbit pAb (A11195) at dilution of 1:50 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate buffer (pH 6.0) prior to IHC staining.

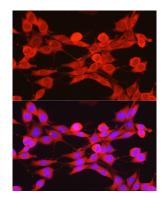


Immunofluorescence analysis of PC-12 cells using FAK Rabbit pAb (A11195) at dilution of 1:50 (40x lens). Secondary antibody: Cy3conjugated Goat anti-Rabbit IgG (H+L) (AS007) at 1:500 dilution. Blue: DAPI for nuclear staining.

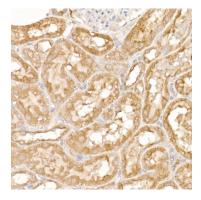


Immunohistochemistry analysis of paraffinembedded Rat kidney using FAK Rabbit pAb (A11195) at dilution of 1:50 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate buffer (pH 6.0) prior to IHC staining.

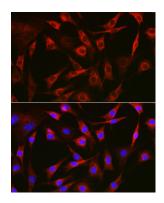


Immunofluorescence analysis of NIH/3T3 cells using FAK Rabbit pAb (A11195) at dilution of 1:100 (40x lens). Secondary antibody: Cy3-conjugated Goat anti-Rabbit IgG (H+L) (AS007) at 1:500 dilution. Blue: DAPI for nuclear staining.

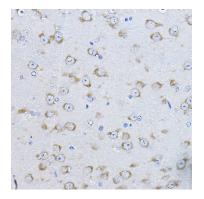

Immunohistochemistry analysis of paraffinembedded Mouse brain using FAK Rabbit pAb (A11195) at dilution of 1:100 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate buffer (pH 6.0) prior to IHC staining.

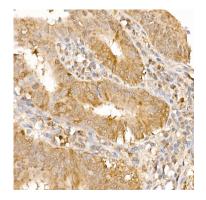


Immunohistochemistry analysis of paraffinembedded Human colon carcinoma using FAK Rabbit pAb (A11195) at dilution of 1:50 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate buffer (pH 6.0) prior to IHC staining.

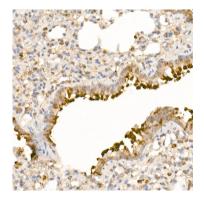


Immunohistochemistry analysis of paraffinembedded Rat lung using FAK Rabbit pAb (A11195) at dilution of 1:50 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate buffer (pH 6.0) prior to IHC staining.




Immunohistochemistry analysis of paraffinembedded Human colon carcinoma using FAK Rabbit pAb (A11195) at dilution of 1:50 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate buffer (pH 6.0) prior to IHC staining.

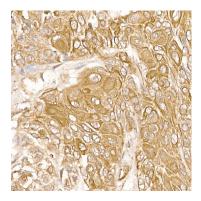


Immunohistochemistry analysis of paraffinembedded Human esophageal cancer using FAK Rabbit pAb (A11195) at dilution of 1:50 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate buffer (pH 6.0) prior to IHC staining.

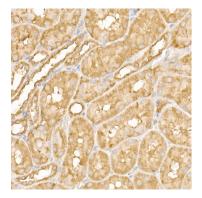


Immunohistochemistry analysis of paraffinembedded Rat kidney using FAK Rabbit pAb (A11195) at dilution of 1:50 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate buffer (pH 6.0) prior to IHC staining.

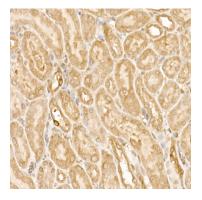


Immunohistochemistry analysis of paraffinembedded Mouse kidney using FAK Rabbit pAb (A11195) at dilution of 1:50 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate buffer (pH 6.0) prior to IHC staining.

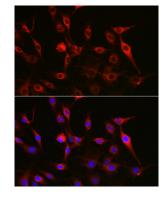


Immunofluorescence analysis of C6 cells using FAK Rabbit pAb (A11195) at dilution of 1:100 (40x lens). Secondary antibody: Cy3conjugated Goat anti-Rabbit IgG (H+L) (AS007) at 1:500 dilution. Blue: DAPI for nuclear staining.

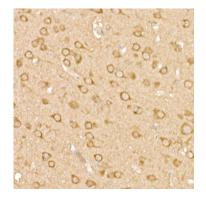


Immunohistochemistry analysis of paraffinembedded Rat brain using FAK Rabbit pAb (A11195) at dilution of 1:50 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate buffer (pH 6.0) prior to IHC staining.

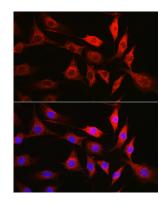


Immunofluorescence analysis of NIH/3T3 cells using FAK Rabbit pAb (A11195) at dilution of 1:100 (40x lens). Secondary antibody: Cy3-conjugated Goat anti-Rabbit IgG (H+L) (AS007) at 1:500 dilution. Blue: DAPI for nuclear staining.

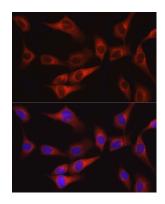

Immunofluorescence analysis of NIH/3T3 cells using FAK Rabbit pAb (A11195) at dilution of 1:100 (40x lens). Secondary antibody: Cy3-conjugated Goat anti-Rabbit IgG (H+L) (AS007) at 1:500 dilution. Blue: DAPI for nuclear staining.

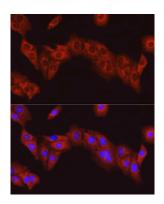

Immunofluorescence analysis of PC-12 cells using FAK Rabbit pAb (A11195) at dilution of 1:100 (40x lens). Secondary antibody: Cy3conjugated Goat anti-Rabbit IgG (H+L) (AS007) at 1:500 dilution. Blue: DAPI for nuclear staining.

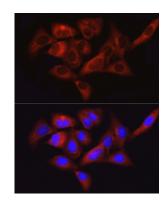

Immunofluorescence analysis of U2OS cells using FAK Rabbit pAb (A11195) at dilution of 1:100 (40x lens). Secondary antibody: Cy3conjugated Goat anti-Rabbit IgG (H+L) (AS007) at 1:500 dilution. Blue: DAPI for nuclear staining.