ACTN3 Rabbit pAb

Catalog No.: A12797



Basic Information

Observed MW

103kDa

Calculated MW

103kDa

Category

Primary antibody

Applications

ELISA,WB,IHC-P

Cross-Reactivity

Human, Mouse, Rat

Background

This gene encodes a member of the alpha-actin binding protein gene family. The encoded protein is primarily expressed in skeletal muscle and functions as a structural component of sarcomeric Z line. This protein is involved in crosslinking actin containing thin filaments. An allelic polymorphism in this gene results in both coding and non-coding variants; the reference genome represents the coding allele. The non-functional allele of this gene is associated with elite athlete status.

Recommended Dilutions

WB 1:500 - 1:2000

IHC-P 1:100 - 1:500

Immunogen Information

Gene ID89

Swiss Prot
Q08043

Immunogen

Recombinant fusion protein containing a sequence corresponding to amino acids 580-630 of human ACTN3 (NP_001095.2).

Synonyms

ACTN3D; ACTN3

Contact

a	400-999-6126
×	cn.market@abclonal.com.cn
$\overline{\Box}$	www abclonal com cn


Product Information

SourceIsotypePurificationRabbitIgGAffinity purification

Storage

Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS with 0.01% thimerosal,50% glycerol,pH7.3.

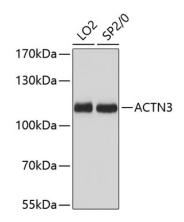

Western blot analysis of various lysates using ACTN3 Rabbit pAb (A12797) at 1:3000 dilution. Secondary antibody: HRP Goat Anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.

Lysates/proteins: 25µg per lane.


Blocking buffer: 3% nonfat dry milk in TBST.

Detection: ECL Basic Kit (RM00020).

Exposure time: 10s.



Immunohistochemistry analysis of paraffinembedded mouse skeletal muscle using ACTN3 Rabbit pAb (A12797) at dilution of 1:500 (40x lens).Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.

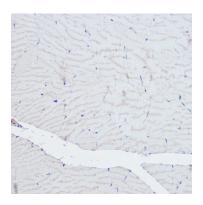

Immunohistochemistry analysis of paraffinembedded rat skeletal muscle using ACTN3 Rabbit pAb (A12797) at dilution of 1:500 (40x lens).Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.