Tyrosinase Rabbit pAb

Catalog No.: A13608 1 Publications



Basic Information

Observed MW

80kDa

Calculated MW

60kDa

Category

Primary antibody

Applications

WB,IHC-P,ELISA

Cross-Reactivity

Human, Rat

Background

The enzyme encoded by this gene catalyzes the first 2 steps, and at least 1 subsequent step, in the conversion of tyrosine to melanin. The enzyme has both tyrosine hydroxylase and dopa oxidase catalytic activities, and requires copper for function. Mutations in this gene result in oculocutaneous albinism, and nonpathologic polymorphisms result in skin pigmentation variation. The human genome contains a pseudogene similar to the 3' half of this gene.

Recommended Dilutions

WB 1:500 - 1:2000

IHC-P 1:1000 - 1:5000

ELISA Recommended starting concentration is 1 μg/mL.

Please optimize the concentration based on your specific assay requirements.

Immunogen Information

Gene ID7299

Swiss Prot
P14679

Immunogen

Recombinant protein (or fragment). This information is considered to be commercially sensitive.

Synonyms

ATN; CMM8; OCA1; OCA1A; OCAIA; SHEP3; Tyrosinase

Contact

6		400-999-6126
\bowtie		cn.market@abclonal.com.cn
•	T	www.abclonal.com.cn


Product Information

SourceIsotypePurificationRabbitIgGAffinity purification

Storage

Store at -20 $^{\circ}\text{C}.$ Avoid freeze / thaw cycles.

Buffer: Buffer: PBS containing 50% glycerol, preserved with proclin300 or sodium azide (as specified on the Certificate of Analysis), pH 7.3.

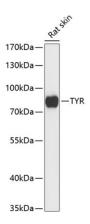


Western blot analysis of lysates from rat skin, using Tyrosinase Rabbit pAb (A13608) at 1:3000 dilution. Secondary antibody: HRP-conjugated Goat anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.

Lysates/proteins: 25µg per lane.

Blocking buffer: 3% nonfat dry milk in TBST.

Detection: ECL Basic Kit (RM00020).

Exposure time: 10s.

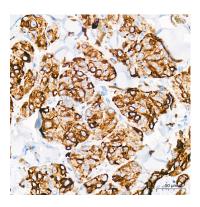

Immunohistochemistry analysis of paraffinembedded Human malignant melanoma tissue using Tyrosinase Rabbit pAb (A13608) at a dilution of 1:2000 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate buffer (pH 6.0) prior to IHC staining.