CCDC124 Rabbit pAb

Catalog No.: A14436



Basic Information

Observed MW

26kDa

Calculated MW

26kDa

Category

Primary antibody

Applications

WB,IHC-P,ELISA

Cross-Reactivity

Human, Mouse, Rat

Background

Enables RNA binding activity. Predicted to be involved in cell division. Located in cytosol and plasma membrane.

Recommended Dilutions

WB 1:500 - 1:2000

IHC-P 1:50 - 1:200

ELISA Recommended starting concentration is 1 μ g/mL. Please optimize the

concentration based on your specific assay requirements.

Immunogen Information

Gene IDSwiss Prot
115098
Q96CT7

Immunogen

Recombinant protein (or fragment). This information is considered to be commercially sensitive.

Synonyms

Lso2; oxs1; CCDC124

Contact

<u>a</u>		400-999-6126
\bowtie		cn.market@abclonal.com.cn
$\overline{\triangle}$	ī	www.ahclonal.com.cn


Product Information

SourceIsotypePurificationRabbitIgGAffinity purification

Storage

Store at -20 $^{\circ}\text{C}.$ Avoid freeze / thaw cycles.

Buffer: PBS with 0.01% thimerosal,50% glycerol,pH7.3.

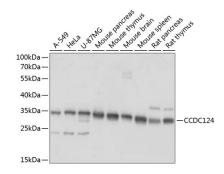


Western blot analysis of various lysates using CCDC124 Rabbit pAb (A14436) at 1:1000 dilution. Secondary antibody: HRP-conjugated Goat anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution. Lysates/proteins: $25\mu g$ per lane.


Blocking buffer: 3% nonfat dry milk in TBST.

Detection: ECL Basic Kit (RM00020).


Exposure time: 1s.

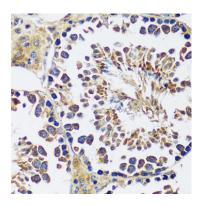

Immunohistochemistry analysis of paraffinembedded Mouse testis using CCDC124 Rabbit pAb (A14436) at dilution of 1:100 (40x lens). Microwave antigen retrieval performed with 0.01M PBS Buffer (pH 7.2) prior to IHC staining.

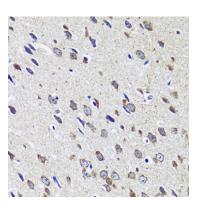

Immunohistochemistry analysis of paraffinembedded Mouse brain using CCDC124 Rabbit pAb (A14436) at dilution of 1:100 (40x lens). Microwave antigen retrieval performed with 0.01M PBS Buffer (pH 7.2) prior to IHC staining.

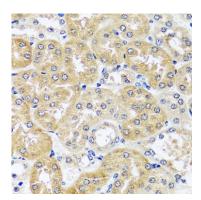

Immunohistochemistry analysis of paraffinembedded Mouse kidney using CCDC124 Rabbit pAb (A14436) at dilution of 1:100 (40x lens). Microwave antigen retrieval performed with 0.01M PBS Buffer (pH 7.2) prior to IHC staining.