# TRPV6 Rabbit pAb

Catalog No.: A16128 4 Publications



### **Basic Information**

#### **Observed MW**

70kDa

#### **Calculated MW**

87kDa

#### Category

Primary antibody

### **Applications**

ELISA,WB,IHC-P

#### **Cross-Reactivity**

Human, Mouse, Rat

### **Background**

This gene encodes a member of a family of multipass membrane proteins that functions as calcium channels. The encoded protein contains N-terminal ankyrin repeats, which are required for channel assembly and regulation. Translation initiation for this protein occurs at a non-AUG start codon that is decoded as methionine. This gene is situated next to a closely related gene for transient receptor potential cation channel subfamily V member 5 (TRPV5). This locus has experienced positive selection in non-African populations, resulting in several non-synonymous codon differences among individuals of different genetic backgrounds.

### **Recommended Dilutions**

**WB** 1:500 - 1:1000

**IHC-P** 1:50 - 1:200

## **Immunogen Information**

**Gene ID**55503

Swiss Prot
Q9H1D0

#### **Immunogen**

A synthetic peptide corresponding to a sequence within amino acids 1-100 of human TRPV6 (NP $\_$ 061116.5).

### **Synonyms**

CAT1; CATL; ZFAB; ECAC2; ABP/ZF; HRPTTN; LP6728; HSA277909; TRPV6

### **Contact**

6	400-999-6126
$\bowtie$	cn.market@abclonal.com.cn
•	www.abclonal.com.cn


### **Product Information**

SourceIsotypePurificationRabbitIgGAffinity purification

#### **Storage**

Store at -20  $^{\circ}\text{C}.$  Avoid freeze / thaw cycles.

Buffer: PBS with 0.01% thimerosal,50% glycerol,pH7.3.

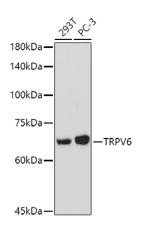


Western blot analysis of various lysates using TRPV6 Rabbit pAb (A16128) at 1:1000 dilution. Secondary antibody: HRP Goat Anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.

Lysates/proteins: 25µg per lane.

Blocking buffer: 3% nonfat dry milk in TBST.

Detection: ECL Basic Kit (RM00020).

Exposure time: 90s.

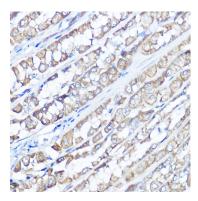

Immunohistochemistry analysis of paraffinembedded rat stomach using TRPV6 Rabbit pAb (A16128) at dilution of 1:100 (40x lens).Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.