SPARC Rabbit pAb

Catalog No.: A1615



Basic Information

Observed MW

35kDa

Calculated MW

35kDa

Category

Primary antibody

Applications

IHC-P,ELISA

Cross-Reactivity

Human, Mouse, Rat

Background

This gene encodes a cysteine-rich acidic matrix-associated protein. The encoded protein is required for the collagen in bone to become calcified but is also involved in extracellular matrix synthesis and promotion of changes to cell shape. The gene product has been associated with tumor suppression but has also been correlated with metastasis based on changes to cell shape which can promote tumor cell invasion. Three transcript variants encoding different isoforms have been found for this gene.

Recommended Dilutions

IHC-P

1:50 - 1:200

ELISA

Recommended starting concentration is 1 µg/mL.
Please optimize the concentration based on your specific assay requirements.

Immunogen Information

Gene ID 6678 Swiss Prot P09486

ImmunogenRecombinant protein (or fragment). This information is considered to be commercially sensitive.

Synonyms

ON; ONT; OI17; BM-40; SPARC

Contact

a	400-999-6126
\sim	cn.market@abclonal.com.cn
⊙	www.abclonal.com.cn

Product Information


SourceIsotypePurificationRabbitIgGAffinity purification

Storage


Store at -20 $^{\circ}\text{C}.$ Avoid freeze / thaw cycles.

Buffer: PBS with 0.09% Sodium azide,50% glycerol,pH7.3.


Validation Data

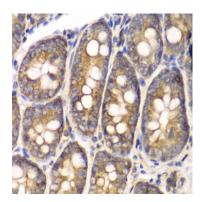

Immunohistochemistry analysis of paraffinembedded Rat Intestine using SPARC Rabbit pAb (A1615) at dilution of 1:100 (40x lens). Microwave antigen retrieval performed with 0.01M PBS Buffer (pH 7.2) prior to IHC staining.

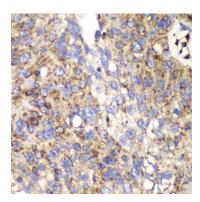

Immunohistochemistry analysis of paraffinembedded Human lung cancer using SPARC Rabbit pAb (A1615) at dilution of 1:100 (40x lens). Microwave antigen retrieval performed with 0.01M PBS Buffer (pH 7.2) prior to IHC staining.

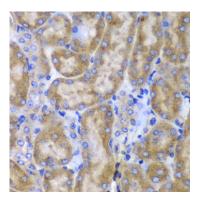

Immunohistochemistry analysis of paraffinembedded Mouse kidney using SPARC Rabbit pAb (A1615) at dilution of 1:100 (40x lens). Microwave antigen retrieval performed with 0.01M PBS Buffer (pH 7.2) prior to IHC staining.