CAMP Rabbit pAb

Catalog No.: A17298



Basic Information

Observed MW

Refer to figures

Calculated MW

19kDa

Category

Primary antibody

Applications

ELISA,IHC-P,IF/ICC

Cross-Reactivity

Human, Rat

Background

This gene encodes a member of an antimicrobial peptide family, characterized by a highly conserved N-terminal signal peptide containing a cathelin domain and a structurally variable cationic antimicrobial peptide, which is produced by extracellular proteolysis from the C-terminus. The protein plays an important role in innate immunity defense against viruses. In addition to its antibacterial, antifungal, and antiviral activities, the encoded protein functions in cell chemotaxis, immune mediator induction, and inflammatory response regulation.

Recommended Dilutions

IHC-P	1:50 - 1:200
IF/ICC	1:50 - 1:200

Immunogen Information

Gene ID820

Swiss Prot
P49913

Immunogen

A synthetic peptide corresponding to a sequence within amino acids 100-170 of human CAMP (NP_004336.4).

Synonyms

LL37; CAP18; CRAMP; HSD26; CAP-18; FALL39; FALL-39; CAMP

Contact

6	400-999-6126
\bowtie	cn.market@abclonal.com.cn
\odot	www.abclonal.com.cn


Product Information

SourceIsotypePurificationRabbitIgGAffinity purification


Storage

Store at -20°C. Avoid freeze / thaw cycles.


Buffer: PBS with 0.01% thimerosal,50% glycerol,pH7.3.

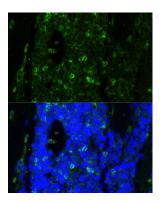

Immunofluorescence analysis of paraffinembedded rat bone marrow using CAMP Rabbit pAb (A17298) at dilution of 1:100. Secondary antibody: Cy3 Goat Anti-Rabbit IgG (H+L) (AS007) at 1:500 dilution. Blue: DAPI for nuclear staining.

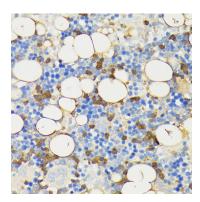

Immunohistochemistry analysis of CAMP in paraffin-embedded rat bone marrow using CAMP Rabbit pAb (A17298) at dilution of 1:200 (40x lens).Perform microwave antigen retrieval with 10 mM PBS buffer pH 7.2 before commencing with IHC staining protocol.

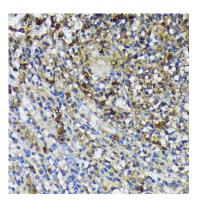

Immunohistochemistry analysis of CAMP in paraffin-embedded human appendix using CAMP Rabbit pAb (A17298) at dilution of 1:200 (40x lens).Perform microwave antigen retrieval with 10 mM PBS buffer pH 7.2 before commencing with IHC staining protocol.