# Ferritin Light Chain Rabbit pAb

Catalog No.: A1768 3 Publications



# **Basic Information**

#### **Observed MW**

20kDa

### **Calculated MW**

20kDa

### Category

Primary antibody

## **Applications**

ELISA, WB, IHC-P

#### **Cross-Reactivity**

Human, Mouse, Rat

# **Background**

This gene encodes the light subunit of the ferritin protein. Ferritin is the major intracellular iron storage protein in prokaryotes and eukaryotes. It is composed of 24 subunits of the heavy and light ferritin chains. Variation in ferritin subunit composition may affect the rates of iron uptake and release in different tissues. A major function of ferritin is the storage of iron in a soluble and nontoxic state. Defects in this light chain ferritin gene are associated with several neurodegenerative diseases and hyperferritinemia-cataract syndrome. This gene has multiple pseudogenes.

## **Recommended Dilutions**

**WB** 1:500 - 1:1000

**IHC-P** 1:50 - 1:200

# Immunogen Information

**Gene ID**2512

Swiss Prot
P02792

#### Immunogen

Recombinant fusion protein containing a sequence corresponding to amino acids 1-175 of human Ferritin Light Chain (NP\_000137.2).

## **Synonyms**

LFTD; NBIA3; Ferritin Light Chain

## **Contact**

<b>a</b>	400-999-6126
×	cn.market@abclonal.com.cn
	www.abclonal.com.cn


## **Product Information**

SourceIsotypePurificationRabbitIgGAffinity purification

#### Storage

Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS with 0.05% proclin300,50% glycerol,pH7.3.

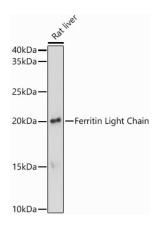

Western blot analysis of lysates from Rat liver, using Ferritin Light Chain Rabbit pAb (A1768) at 1:1000 dilution


Secondary antibody: HRP Goat Anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.

Lysates/proteins: 25µg per lane.

Blocking buffer: 3% nonfat dry milk in TBST.

Detection: ECL Basic Kit (RM00020).

Exposure time: 180s.

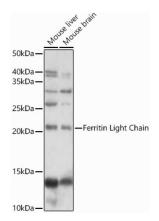


Western blot analysis of various lysates using Ferritin Light Chain Rabbit pAb (A1768) at 1:1000 dilution. Secondary antibody: HRP Goat Anti-Rabbit  $\lg G$  (H+L) (AS014) at 1:10000 dilution.

Lysates/proteins: 25µg per lane.


Blocking buffer: 3% nonfat dry milk in TBST.

Detection: ECL Basic Kit (RM00020).


Exposure time: 30s.

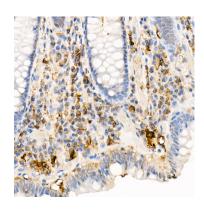

Immunohistochemistry analysis of paraffinembedded human colon using Ferritin Light Chain Rabbit pAb (A1768) at dilution of 1:100 (40x lens).Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.

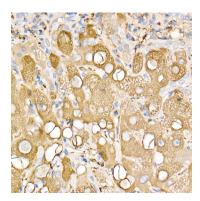


Immunohistochemistry analysis of paraffinembedded human liver using Ferritin Light Chain Rabbit pAb (A1768) at dilution of 1:100 (40x lens).Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.

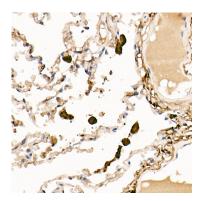

Immunohistochemistry analysis of paraffinembedded human lung using Ferritin Light Chain Rabbit pAb (A1768) at dilution of 1:100 (40x lens).Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol


## **Validation Data**

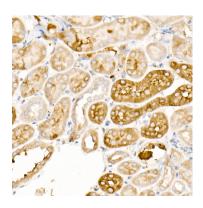


Immunohistochemistry analysis of paraffinembedded mouse kidney using Ferritin Light Chain Rabbit pAb (A1768) at dilution of 1:100 (40x lens).Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol



Immunohistochemistry analysis of paraffinembedded mouse liver using Ferritin Light Chain Rabbit pAb (A1768) at dilution of 1:100 (40x lens).Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.

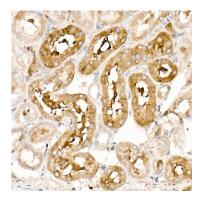

Immunohistochemistry analysis of paraffinembedded rat kidney using Ferritin Light Chain Rabbit pAb (A1768) at dilution of 1:100 (40x lens).Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol

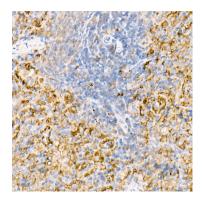

Immunohistochemistry analysis of paraffinembedded rat spleen using Ferritin Light Chain Rabbit pAb (A1768) at dilution of 1:100 (40x lens).Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.