SENP2 Rabbit pAb

Catalog No.: A17994



Basic Information

Observed MW

68kDa

Calculated MW

68kDa

Category

Primary antibody

Applications

ELISA,WB,IHC-P

Cross-Reactivity

Human, Mouse, Rat

Background

SUMO1 (UBL1; MIM 601912) is a small ubiquitin-like protein that can be covalently conjugated to other proteins. SENP2 is one of a group of enzymes that process newly synthesized SUMO1 into the conjugatable form and catalyze the deconjugation of SUMO1-containing species.

Recommended Dilutions

WB 1:500 - 1:2000

IHC-P 1:50 - 1:200

Immunogen Information

Gene IDSwiss Prot
59343
Q9HC62

Immunogen

Recombinant fusion protein containing a sequence corresponding to amino acids 366-589 of human SENP2 (NP_067640.2).

Synonyms

AXAM2; SMT3IP2; SENP2

Contact

a		400-999-6126
\bowtie		cn.market@abclonal.com.cn
\overline{a}	ı	www.ahclonal.com.cn


Product Information

SourceIsotypePurificationRabbitIgGAffinity purification

Storage

Store at -20 $^{\circ}\text{C}.$ Avoid freeze / thaw cycles.

Buffer: PBS with 0.01% thimerosal,50% glycerol,pH7.3.

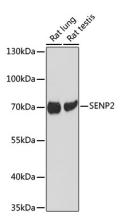


Western blot analysis of various lysates, using SENP2 Rabbit pAb (A17994) at 1:1000 dilution. Secondary antibody: HRP Goat Anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.


Lysates/proteins: 25µg per lane.

Blocking buffer: 3% nonfat dry milk in TBST. Detection: ECL Basic Kit (RM00020).


Exposure time: 90s.

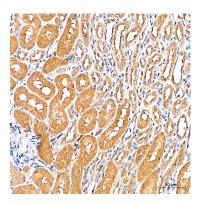

Immunohistochemistry analysis of paraffinembedded Mouse kidney tissue using SENP2 Rabbit pAb (A17994) at a dilution of 1:100 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate Buffer(pH 6.0) prior to IHC staining.

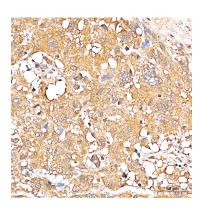

Immunohistochemistry analysis of paraffinembedded Human lung cancer tissue using SENP2 Rabbit pAb (A17994) at a dilution of 1:100 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate Buffer(pH 6.0) prior to IHC staining.

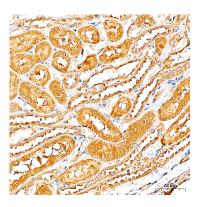

Immunohistochemistry analysis of paraffinembedded Rat kidney tissue using SENP2 Rabbit pAb (A17994) at a dilution of 1:100 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate Buffer(pH 6.0) prior to IHC staining.