# PTGIR Rabbit pAb

Catalog No.: A1849



## **Basic Information**

#### **Observed MW**

41kDa

### **Calculated MW**

41kDa

### Category

Primary antibody

### **Applications**

ELISA,WB,IHC-P

#### **Cross-Reactivity**

Human, Mouse, Rat

# **Background**

The protein encoded by this gene is a member of the G-protein coupled receptor family 1 and has been shown to be a receptor for prostacyclin. Prostacyclin, the major product of cyclooxygenase in macrovascular endothelium, elicits a potent vasodilation and inhibition of platelet aggregation through binding to this receptor.

# **Recommended Dilutions**

**WB** 1:500 - 1:2000

**IHC-P** 1:50 - 1:200

# **Immunogen Information**

**Gene ID**5739

Swiss Prot
P43119

#### **Immunogen**

Recombinant fusion protein containing a sequence corresponding to amino acids 207-386 of human PTGIR (NP\_000951.1).

### **Synonyms**

IP; PRIPR; PTGIR

### **Contact**

<u>a</u>		400-999-6126
$\bowtie$		cn.market@abclonal.com.cn
$\overline{a}$	ı	www.ahclonal.com.cn


### **Product Information**

SourceIsotypePurificationRabbitIgGAffinity purification

### **Storage**

Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS with 0.02% sodium azide,50% glycerol,pH7.3.

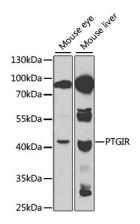



Western blot analysis of various lysates using PTGIR Rabbit pAb (A1849) at 1:1000 dilution. Secondary antibody: HRP Goat Anti-Rabbit  $\lg G$  (H+L) (AS014) at 1:10000 dilution.

Lysates/proteins: 25µg per lane.

Blocking buffer: 3% nonfat dry milk in TBST. Detection: ECL Enhanced Kit (RM00021).

Exposure time: 90s.

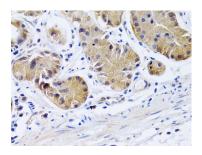

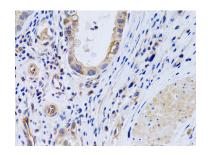

Immunohistochemistry analysis of paraffinembedded human stomach using PTGIR Rabbit pAb (A1849) at dilution of 1:200 (40x lens).Perform microwave antigen retrieval with 10 mM PBS buffer pH 7.2 before commencing with IHC staining protocol.

Immunohistochemistry analysis of paraffinembedded human gastric cancer using PTGIR Rabbit pAb (A1849) at dilution of 1:200 (40x lens).Perform microwave antigen retrieval with 10 mM PBS buffer pH 7.2 before commencing with IHC staining protocol.