WWC2 Rabbit pAb

Catalog No.: A18668



Basic Information

Observed MW

134kDa

Calculated MW

134kDa

Category

Primary antibody

Applications

ELISA,WB,IHC-P

Cross-Reactivity

Human, Mouse, Rat

Background

This gene encodes a member of the WW-and-C2-domain-containing family of proteins. Members of this family have two N-terminal WW domains that mediate binding to target proteins harboring L/PPxY motifs, an internal C2 domain for membrane association, and C-terminal alpha protein kinase C binding sites and class III PDZ domain-interaction motifs. Proteins of this family are able to form homo- and heterodimers and to modulate hippo pathway signaling.

Recommended Dilutions

WB 1:500 - 1:2000

IHC-P 1:50 - 1:200

Immunogen Information

Gene ID80014

Swiss Prot
Q6AWC2

Immunogen

Recombinant fusion protein containing a sequence corresponding to amino acids 160-430 of human WWC2 (NP_079225.5).

Synonyms

BOMB; WWC2

Contact

<u>a</u>		400-999-6126
\bowtie		cn.market@abclonal.com.cn
$\overline{\triangle}$	ī	www.ahclonal.com.cn


Product Information

SourceIsotypePurificationRabbitIgGAffinity purification

Storage

Store at -20 $^{\circ}\text{C}.$ Avoid freeze / thaw cycles.

Buffer: PBS with 0.01% thimerosal,50% glycerol,pH7.3.

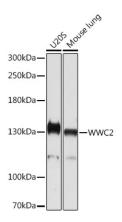


Western blot analysis of various lysates using WWC2 Rabbit pAb (A18668) at 1:1000 dilution. Secondary antibody: HRP Goat Anti-Rabbit lgG (H+L) (AS014) at 1:10000 dilution. Lysates/proteins: $25\mu g$ per lane.


Blocking buffer: 3% nonfat dry milk in TBST.

Detection: ECL Basic Kit (RM00020).

Exposure time: 5s.

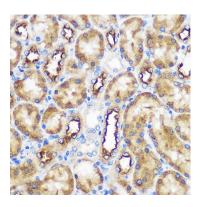

Immunohistochemistry analysis of WWC2 in paraffin-embedded Rat kidney using WWC2 Rabbit pAb (A18668) at dilution of 1:100 (40x lens).Perform microwave antigen retrieval with 10 mM PBS buffer pH 7.2 before commencing with IHC staining protocol.

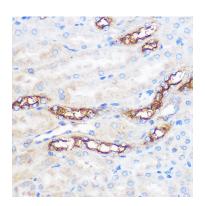

Immunohistochemistry analysis of WWC2 in paraffin-embedded Mouse kidney using WWC2 Rabbit pAb (A18668) at dilution of 1:100 (40x lens).Perform microwave antigen retrieval with 10 mM PBS buffer pH 7.2 before commencing with IHC staining protocol.