# CD34 Rabbit mAb

Catalog No.: A19015 Recombinant 10 Publications



## **Basic Information**

#### **Observed MW**

120kDa

### **Calculated MW**

41kDa

### Category

Primary antibody

### **Applications**

WB,IHC-P,ELISA

### **Cross-Reactivity**

Human, Mouse, Rat

#### CloneNo number

ARC0219

## **Background**

The protein encoded by this gene may play a role in the attachment of stem cells to the bone marrow extracellular matrix or to stromal cells. This single-pass membrane protein is highly glycosylated and phosphorylated by protein kinase C. Two transcript variants encoding different isoforms have been found for this gene.

## **Recommended Dilutions**

WB 1:500 - 1:1000

**IHC-P** 1:100 - 1:400

**ELISA** Recommended starting concentration is 1 μg/mL.

Please optimize the concentration based on your specific assay requirements.

## **Immunogen Information**

**Gene ID**947

Swiss Prot
P28906

### **Immunogen**

Synthetic peptide. This information is considered to be commercially sensitive.

### **Synonyms**

CD34; CD34 molecule; GIG3; MORT1

## **Contact**

<b>a</b>	400-999-6126
$\sim$	cn.market@abclonal.com.cn
$\overline{\Box}$	www.ahclonal.com.cn


## **Product Information**

SourceIsotypePurificationRabbitIgGAffinity purification

#### Storage

Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS with 0.02% sodium azide, 0.05% BSA, 50% glycerol, pH7.3.

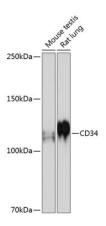


Western blot analysis of various lysates using CD34 Rabbit mAb (A19015) at 1:500 dilution. Secondary antibody: HRP-conjugated Goat anti-Rabbit lgG (H+L) (AS014) at 1:10000 dilution. Lysates/proteins:  $25\mu g$  per lane.


Blocking buffer: 3% nonfat dry milk in TBST.

Detection: ECL Basic Kit (RM00020).


Exposure time: 3min.

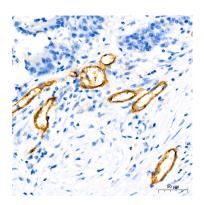


Immunohistochemistry analysis of paraffinembedded Human colon carcinoma using CD34 Rabbit mAb (A19015) at dilution of 1:100 (40x lens). Microwave antigen retrieval performed with 0.01M Tris/EDTA Buffer (pH 9.0) prior to IHC staining.

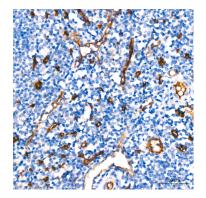


Immunohistochemistry analysis of paraffinembedded Human tonsil using CD34 Rabbit mAb (A19015) at dilution of 1:100 (40x lens). Microwave antigen retrieval performed with 0.01M Tris/EDTA Buffer (pH 9.0) prior to IHC staining.

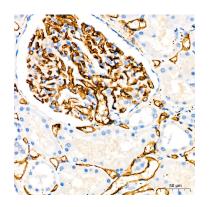


Immunohistochemistry analysis of paraffinembedded Human kidney using CD34 Rabbit mAb (A19015) at dilution of 1:100 (40x lens). Microwave antigen retrieval performed with 0.01M Tris/EDTA Buffer (pH 9.0) prior to IHC staining.

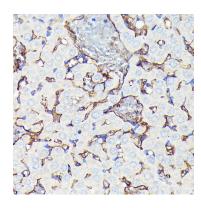

Immunohistochemistry analysis of paraffinembedded Rat ovary tissue using CD34 Rabbit mAb (A19015) at a dilution of 1:100 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer (pH 9.0) prior to IHC staining.

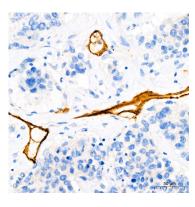


Immunohistochemistry analysis of paraffinembedded Human lung cancer using CD34 Rabbit mAb (A19015) at dilution of 1:100 (40x lens). Microwave antigen retrieval performed with 0.01M Tris/EDTA Buffer (pH 9.0) prior to IHC staining.

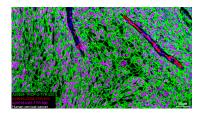

The multiplex IHC analysis on paraffinembedded Human cervical cancer tissue using the following specific primary antibodies and tyramide signal amplification (TSA) reagents (RK05903): TROP-2 Rabbit mAb (A20824, 1:8000) with TSA-TYR-520 (Green), CD34 Rabbit mAb (A19015, 1:100) with TSA-TYR-570 (Red), and Ki67 Rabbit mAb (A20018, 1:500) with TSA-TYR-690 (Magenta). DAPI (Blue) was used for nuclear staining. Prior to multiplex IHC staining, high-pressure antigen retrieval was performed using 0.01M citrate buffer at pH 6.0. The analysis was completed using a 20x

objective lens.

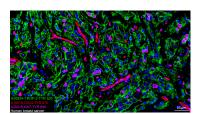

The multiplex IHC analysis on paraffinembedded Human breast cancer tissue using the following specific primary antibodies and tyramide signal amplification (TSA) reagents (RK05903): TROP-2 Rabbit mAb (A20824, 1:8000) with TSA-TYR-520 (Green), CD34 Rabbit mAb (A19015, 1:100) with TSA-TYR-570 (Red), and Ki67 Rabbit mAb (A20018, 1:500) with TSA-TYR-690 (Magenta). DAPI (Blue) was used for nuclear staining. Prior to multiplex IHC staining, highpressure antigen retrieval was performed using 0.01M citrate buffer at pH 6.0. The analysis was completed using a 20x objective lens.