KAP1/TRIM28 Rabbit mAb

Catalog No.: A19568 Recombinant 2 Publications



Basic Information

Observed MW

110kDa

Calculated MW

89kDa

Category

Primary antibody

Applications

WB,IHC-P,IP,ELISA

Cross-Reactivity

Human, Mouse, Rat

CloneNo number

ARC0047

Background

The protein encoded by this gene mediates transcriptional control by interaction with the Kruppel-associated box repression domain found in many transcription factors. The protein localizes to the nucleus and is thought to associate with specific chromatin regions. The protein is a member of the tripartite motif family. This tripartite motif includes three zincbinding domains, a RING, a B-box type 1 and a B-box type 2, and a coiled-coil region.

Recommended Dilutions

WB 1:1000 - 1:6000

IHC-P 1:200 - 1:2000

IP 0.5μg-4μg antibody for

200μg-400μg extracts of

whole cells

ELISA Recommended starting

concentration is 1 µg/mL.

Please optimize the
concentration based on
your specific assay
requirements.

Immunogen Information

Gene ID Swiss Prot 10155 Q13263

Immunogen

A synthetic peptide corresponding to a sequence within amino acids 250-350 of human KAP1/TRIM28 (Q13263).

Synonyms

KAP1; TF1B; RNF96; TIF1B; PPP1R157; TIF1beta; KAP1/TRIM28

Contact

| a | 400-999-6126 |
|---------------------------|---------------------------|
| \times | cn.market@abclonal.com.cn |
| $\overline{\mathfrak{S}}$ | www.abclonal.com.cn |


Product Information

SourceIsotypePurificationRabbitIgGAffinity purification


Storage

Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS with 0.02% sodium azide, 0.05% BSA, 50% glycerol, pH7.3.

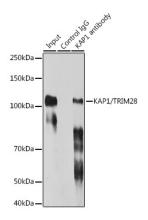

Immunoprecipitation analysis of 300 μg extracts of HeLa cells using 3 μg KAP1/TRIM28 antibody (A19568). Western blot was performed from the immunoprecipitate using KAP1/TRIM28 antibody (A19568) at a dilution of 1:1000.

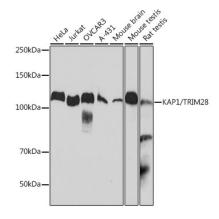


Western blot analysis of various lysates using KAP1/TRIM28 Rabbit mAb (A19568) at 1:1000 dilution. Secondary antibody: HRP-conjugated Goat anti-Rabbit lgG (H+L) (AS014) at 1:10000 dilution.

Lysates/proteins: 25µg per lane.

Blocking buffer: 3% nonfat dry milk in TBST. Detection: ECL Basic Kit (RM00020).

Exposure time: 10s.

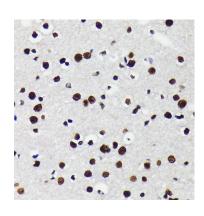


Immunohistochemistry analysis of paraffinembedded Mouse brain using KAP1/TRIM28 Rabbit mAb (A19568) at dilution of 1:100 (40x lens). Microwave antigen retrieval performed with 0.01M PBS Buffer (pH 7.2) prior to IHC staining.



Immunohistochemistry analysis of paraffinembedded Human colon carcinoma using KAP1/TRIM28 Rabbit mAb (A19568) at dilution of 1:200 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate Bufferr (pH 6.0) prior to IHC staining.




Immunohistochemistry analysis of paraffinembedded Human colon using KAP1/TRIM28 Rabbit mAb (A19568) at dilution of 1:200 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate Bufferr (pH 6.0) prior to IHC staining.

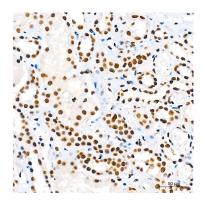


Immunohistochemistry analysis of paraffinembedded Human kidney using KAP1/TRIM28 Rabbit mAb (A19568) at dilution of 1:200 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate Bufferr (pH 6.0) prior to IHC staining.



Immunohistochemistry analysis of paraffinembedded Human spleen using KAP1/TRIM28 Rabbit mAb (A19568) at dilution of 1:200 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate Bufferr (pH 6.0) prior to IHC staining.

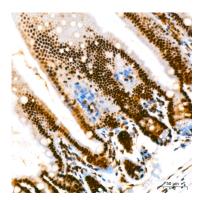

Immunohistochemistry analysis of paraffinembedded Mouse intestin using KAP1/TRIM28 Rabbit mAb (A19568) at dilution of 1:200 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate Bufferr (pH 6.0) prior to IHC staining.

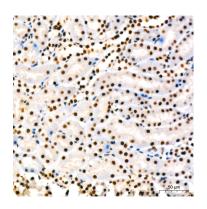

Immunohistochemistry analysis of paraffinembedded Mouse kidney using KAP1/TRIM28 Rabbit mAb (A19568) at dilution of 1:200 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate Bufferr (pH 6.0) prior to IHC staining.



Immunohistochemistry analysis of paraffinembedded Mouse spleen using KAP1/TRIM28 Rabbit mAb (A19568) at dilution of 1:200 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate Bufferr (pH 6.0) prior to IHC staining.