# SOX9 Rabbit mAb

Catalog No.: A19710 Recombinant 28 Publications



# **Basic Information**

#### **Observed MW**

56kDa

### **Calculated MW**

56kDa

### Category

Primary antibody

### **Applications**

ELISA,WB,IHC-P,IF/ICC,IP

#### **Cross-Reactivity**

Human, Mouse, Rat

#### CloneNo number

ARC0190

# **Background**

The protein encoded by this gene recognizes the sequence CCTTGAG along with other members of the HMG-box class DNA-binding proteins. It acts during chondrocyte differentiation and, with steroidogenic factor 1, regulates transcription of the anti-Muellerian hormone (AMH) gene. Deficiencies lead to the skeletal malformation syndrome campomelic dysplasia, frequently with sex reversal.

## **Recommended Dilutions**

**WB** 1:500 - 1:2000

IHC-P 1:100 - 1:500

**IF/ICC** 1:50 - 1:200

**IP** 0.5μg-4μg antibody for

200µg-400µg extracts of

whole cells

# **Immunogen Information**

**Gene ID Swiss Prot**6662
P48436

#### Immunogen

Recombinant fusion protein containing a sequence corresponding to amino acids 169-300 of human SOX9 (P48436).

### **Synonyms**

CMD1; SRA1; CMPD1; SRXX2; SRXY10; SOX9

## **Contact**

<b>a</b>		400-999-6126
$\bowtie$		cn.market@abclonal.com.cn
$\odot$	T	www.abclonal.com.cn


### **Product Information**

SourceIsotypePurificationRabbitIgGAffinity purification


### **Storage**

Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS with 0.02% sodium azide, 0.05% BSA, 50% glycerol, pH7.3.

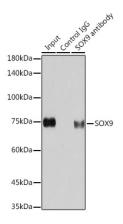

Immunoprecipitation analysis of 200  $\mu g$  extracts of HeLa cells using 3  $\mu g$  SOX9 antibody (A19710). Western blot was performed from the immunoprecipitate using SOX9 antibody (A19710) at a dilution of 1:1000.

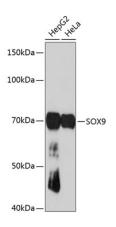


Western blot analysis of various lysates using SOX9 Rabbit mAb (A19710) at 1:1000 dilution. Secondary antibody: HRP Goat Anti-Rabbit  $\lg G$  (H+L) (AS014) at 1:10000 dilution.

Lysates/proteins: 25µg per lane.


Blocking buffer: 3% nonfat dry milk in TBST.

Detection: ECL Basic Kit (RM00020).


Exposure time: 3min.

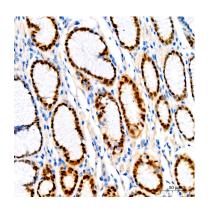


Immunohistochemistry analysis of SOX9 in paraffin-embedded human colon using SOX9 Rabbit mAb (A19710) at dilution of 1:400 (40x lens).Perform high pressure antigen retrieval with 10 mM Tris/EDTA buffer pH 9.0 before commencing with IHC staining protocol.

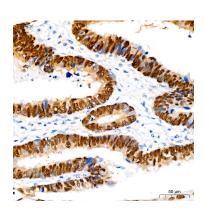


Immunohistochemistry analysis of SOX9 in paraffin-embedded human gastric cancer using SOX9 Rabbit mAb (A19710) at dilution of 1:400 (40x lens).Perform high pressure antigen retrieval with 10 mM Tris/EDTA buffer pH 9.0 before commencing with IHC staining protocol.

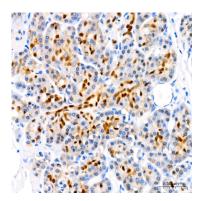


Immunohistochemistry analysis of SOX9 in paraffin-embedded human pancreas using SOX9 Rabbit mAb (A19710) at dilution of 1:400 (40x lens).Perform high pressure antigen retrieval with 10 mM Tris/EDTA buffer pH 9.0 before commencing with IHC staining protocol.

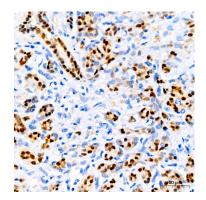


Immunohistochemistry analysis of SOX9 in paraffin-embedded Human pancreatic ductal carcinoma using SOX9 Rabbit mAb (A19710) at dilution of 1:400 (40x lens).Perform high pressure antigen retrieval with 10 mM Tris/EDTA buffer pH 9.0 before commencing with IHC staining protocol.

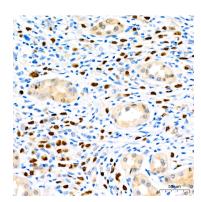


Immunohistochemistry analysis of SOX9 in paraffin-embedded human stomach using SOX9 Rabbit mAb (A19710) at dilution of 1:400 (40x lens).Perform high pressure antigen retrieval with 10 mM Tris/EDTA buffer pH 9.0 before commencing with IHC staining protocol.

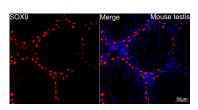

Confocal imaging of paraffin-embedded Mouse testis using SOX9 Rabbit mAb (A19710, dilution 1:200) followed by a further incubation with Cy3 Goat Anti-Rabbit IgG (H+L) (AS007, dilution 1:500) (Red). DAPI was used for nuclear staining (Blue). Objective: 40x.Perform high pressure antigen retrieval with 0.01M citrate buffer (pH 6.0) prior to IF staining.

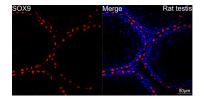

Confocal imaging of paraffin-embedded Rat testis using SOX9 Rabbit mAb (A19710, dilution 1:200) followed by a further incubation with Cy3 Goat Anti-Rabbit IgG (H+L) (AS007, dilution 1:500) (Red). DAPI was used for nuclear staining (Blue). Objective: 40x.Perform high pressure antigen retrieval with 0.01M citrate buffer (pH 6.0) prior to IF staining.

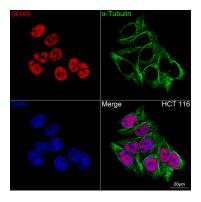

Confocal imaging of HCT 116 cells using SOX9 Rabbit mAb (A19710, dilution 1:200) followed by a further incubation with Cy3 Goat Anti-Rabbit IgG (H+L) (AS007, dilution 1:500) (Red). The cells were counterstained with  $\alpha$ -Tubulin Mouse mAb (AC012, dilution 1:400) followed by incubation with ABflo® 488-conjugated Goat Anti-Mouse IgG (H+L) Ab (AS076, dilution 1:500) (Green). DAPI was used for nuclear staining (Blue). Objective: 100x.