NDUFB8 Rabbit mAb



Catalog No.: A19732 Recombinant 4 Publications

Basic Information

Observed MW

22kDa

Calculated MW

22kDa

Category

Primary antibody

Applications

ELISA,WB,IHC-P,IF/ICC

Cross-Reactivity

Human, Mouse, Rat

CloneNo number

ARC2259

Background

Involved in mitochondrial respiratory chain complex I assembly. Located in endoplasmic reticulum and mitochondrion. Part of mitochondrial respiratory chain complex I. Implicated in nuclear type mitochondrial complex I deficiency 32. Biomarker of Alzheimer's disease and Parkinson's disease.

Recommended Dilutions

WB	1:1000 - 1:3000
IHC-P	1:50 - 1:200
IF/ICC	1.50 - 1.200

Immunogen Information

Gene ID	Swiss Prot
4714	O95169

Immunogen

Recombinant fusion protein containing a sequence corresponding to amino acids 1-150 of human NDUFB8 (095169).

Synonyms

ASHI; CI-ASHI; MC1DN32; NDUFB8

Contact

2		400-999-6126
\bowtie	Τ	cn.market@abclonal.com.cn
•	Т	www.abclonal.com.cn


Product Information

Source	Isotype	Purification
Rabbit	IgG	Affinity purification

Storage

Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS with 0.02% sodium azide, 0.05% BSA, 50% glycerol, pH7.3.

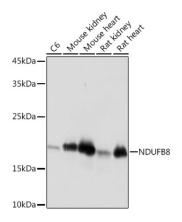


Western blot analysis of various lysates using NDUFB8 Rabbit mAb (A19732) at 1:1000 dilution. Secondary antibody: HRP Goat Anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.

Lysates/proteins: 25µg per lane.

Blocking buffer: 3% nonfat dry milk in TBST.

Detection: ECL Basic Kit (RM00020).

Exposure time: 1s.

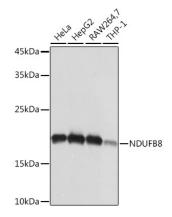


Western blot analysis of various lysates using NDUFB8 Rabbit mAb (A19732) at 1:1000 dilution. Secondary antibody: HRP Goat Anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.


Lysates/proteins: 25µg per lane.


Blocking buffer: 3% nonfat dry milk in TBST.

Detection: ECL Basic Kit (RM00020).

Exposure time: 1s.

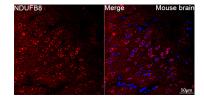

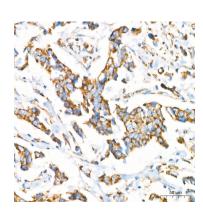


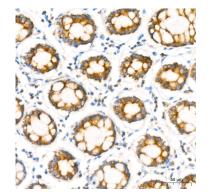


Confocal imaging of paraffin-embedded Mouse brain using NDUFB8 Rabbit mAb (A19732, dilution 1:200) followed by a further incubation with Cy3 Goat Anti-Rabbit IgG (H+L) (AS007, dilution 1:500) (Red). DAPI was used for nuclear staining (Blue). Objective: 40x. Perform microwave antigen retrieval with 0.01 M citrate buffer (pH 6.0) prior to IF staining.

Immunohistochemistry analysis of NDUFB8 in paraffin-embedded human breast cancer using NDUFB8 Rabbit mAb (A19732) at dilution of 1:200 (40x lens).Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.


Immunohistochemistry analysis of NDUFB8 in paraffin-embedded human colon using NDUFB8 Rabbit mAb (A19732) at dilution of 1:200 (40x lens).Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.

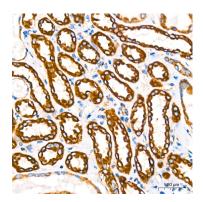


Immunohistochemistry analysis of NDUFB8 in paraffin-embedded human kidney using NDUFB8 Rabbit mAb (A19732) at dilution of 1:200 (40x lens).Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol

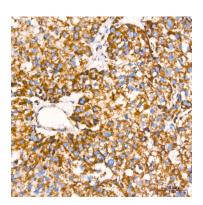


Immunohistochemistry analysis of NDUFB8 in paraffin-embedded human liver cancer using NDUFB8 Rabbit mAb (A19732) at dilution of 1:200 (40x lens).Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.

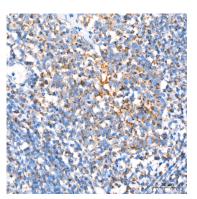

Immunohistochemistry analysis of NDUFB8 in paraffin-embedded human tonsil using NDUFB8 Rabbit mAb (A19732) at dilution of 1:200 (40x lens).Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.

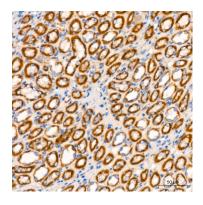

Immunohistochemistry analysis of NDUFB8 in paraffin-embedded rat kidney using NDUFB8 Rabbit mAb (A19732) at dilution of 1:200 (40x lens).Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.

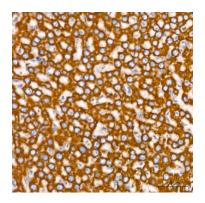

Immunohistochemistry analysis of NDUFB8 in paraffin-embedded rat liver using NDUFB8 Rabbit mAb (A19732) at dilution of 1:200 (40x lens).Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.