Lambda Light chain Rabbit mAb

Catalog No.: A19831 Recombinant



Basic Information

Observed MW

25kDa

Calculated MW

12kDa

Category

Primary antibody

Applications

ELISA,WB,IHC-P

Cross-Reactivity

Human, Mouse, Rat

CloneNo number

ARC2357

Background

Predicted to be involved in immune response. Located in extracellular exosome.

Recommended Dilutions

WB 1:500 - 1:1000

IHC-P 1:50 - 1:200

Immunogen Information

Gene ID28820

Swiss Prot
P01701

Immunogen

Recombinant fusion protein containing a sequence corresponding to amino acids 1-117 of human Lambda Light chain (P01701)

Synonyms

V1-19; IGLV151; Lambda Light chain

Contact

a	400-999-6126
×	cn.market@abclonal.com.cn
	www.abclonal.com.cn


Product Information

SourceIsotypePurificationRabbitIgGAffinity purification

Storage

Store at -20 $^{\circ}\text{C}.$ Avoid freeze / thaw cycles.

Buffer: PBS with 0.02% sodium azide, 0.05% BSA, 50% glycerol, pH7.3.

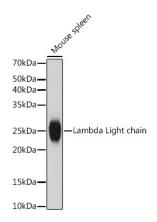

Western blot analysis of lysates from Mouse spleen, using Lambda Light chain Rabbit mAb (A19831) at 1:1000 dilution


Secondary antibody: HRP Goat Anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.

Lysates/proteins: 25µg per lane.

Blocking buffer: 3% nonfat dry milk in TBST.

Detection: ECL Basic Kit (RM00020).

Exposure time: 1s.

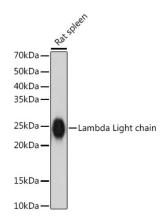

Western blot analysis of lysates from Rat spleen, using Lambda Light chain Rabbit mAb (A19831) at 1:1000 dilution.


Secondary antibody: HRP Goat Anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.

Lysates/proteins: 25µg per lane.


Blocking buffer: 3% nonfat dry milk in TBST.

Detection: ECL Basic Kit (RM00020).

Exposure time: 10s.

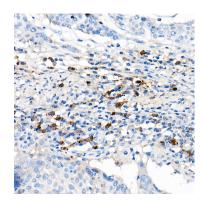

Immunohistochemistry analysis of paraffinembedded human cervix cancer using Lambda Light chain Rabbit mAb (A19831) at dilution of 1:100 (40x lens).Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.

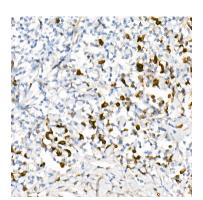

Immunohistochemistry analysis of paraffinembedded human tonsil using Lambda Light chain Rabbit mAb (A19831) at dilution of 1:100 (40x lens).Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.