NDUFB4 Rabbit mAb

Catalog No.: A21115 Recombinant



Basic Information

Observed MW

15kDa

Calculated MW

15kDa

Category

Primary antibody

Applications

ELISA,WB,IHC-P

Cross-Reactivity

Human, Mouse, Rat

CloneNo number

ARC3019

Background

This gene encodes a non-catalytic subunit of the multisubunit NADH:ubiquinone oxidoreductase, the first enzyme complex in the mitochondrial electron transport chain (complex I). Mammalian complex I is composed of 45 different subunits and transfers electrons from NADH to ubiquinone.

Recommended Dilutions

WB 1:500 - 1:1000

IHC-P 1:50 - 1:200

Immunogen Information

Gene ID4710

Swiss Prot

095168

Immunogen

Recombinant fusion protein containing a sequence corresponding to amino acids 1-129 of human NDUFB4 $[NP_004538.2]$.

Synonyms

B15; CI-B15; NDUFB4

Contact

a		400-999-6126
\bowtie		cn.market@abclonal.com.cn
\overline{a}	ı	www.ahclonal.com.cn


Product Information

SourceIsotypePurificationRabbitIgGAffinity purification

Storage

Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS with 0.02% sodium azide, 0.05% BSA, 50% glycerol, pH7.3.

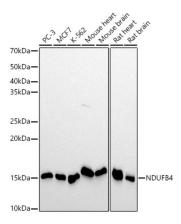


Western blot analysis of various lysates using NDUFB4 Rabbit mAb (A21115) at1:1000 dilution. Secondary antibody: HRP Goat Anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.


Lysates/proteins: 25µg per lane.

Blocking buffer: 3% nonfat dry milk in TBST. Detection: ECL Basic Kit (RM00020).


Exposure time: 30s.

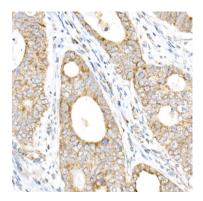

Immunohistochemistry analysis of paraffinembedded human colon carcinoma using NDUFB4 Rabbit mAb (A21115) at dilution of 1:100(40x lens).Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.

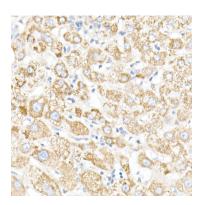

Immunohistochemistry analysis of paraffinembedded human liver using NDUFB4 Rabbit mAb (A21115) at dilution of 1:100(40x lens).Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.

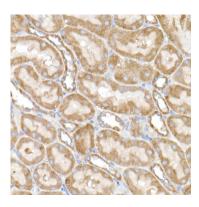

Immunohistochemistry analysis of paraffinembedded rat kidney using NDUFB4 Rabbit mAb (A21115) at dilution of 1:100(40x lens).Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.