# Acetyl-Histone H3-K27 Rabbit mAb

www.abclonal.com

ABclonal

Catalog No.: A22054 Recombinant 1 Publications

## **Basic Information**

### **Observed MW**

17kDa

### **Calculated MW**

16kDa

## Category

Primary antibody

## **Applications**

ELISA, WB, IHC-P, IP, ChIP, ChIPseq,CUT&Tag

## **Cross-Reactivity**

Human, Mouse, Rat, Other (Wide Range Predicted)

## CloneNo number

ARC53902

## **Background**

Histones are basic nuclear proteins that are responsible for the nucleosome structure of the chromosomal fiber in eukaryotes. This structure consists of approximately 146 bp of DNA wrapped around a nucleosome, an octamer composed of pairs of each of the four core histones (H2A, H2B, H3, and H4). The chromatin fiber is further compacted through the interaction of a linker histone, H1, with the DNA between the nucleosomes to form higher order chromatin structures. This gene is intronless and encodes a replication-dependent histone that is a member of the histone H3 family. Transcripts from this gene lack polyA tails; instead, they contain a palindromic termination element. This gene is found in the large histone gene cluster on chromosome 6p22-p21.3.

## **Recommended Dilutions**

**WB** 1:2000 - 1:20000

1:50 - 1:200 **IHC-P** 

0.5μg-4μg antibody for ΙP 200μg-400μg extracts of

whole cells

5µg antibody for ChIP

5μg-10μg of Chromatin

1:20 - 1:100 ChIP-seq

10<sup>5</sup> cells /1 μg **CUT&Tag** 

Contact

<b>a</b>	400-999-6126
$\bowtie$	cn.market@abclonal.com.cn
•	www.abclonal.com.cn

## **Immunogen Information**

**Gene ID Swiss Prot** 8290/8350 Q16695/P68431

#### **Immunogen**

A synthetic acetylated peptide around K4 & K27 of human Histone (NP\_003520.1).

## **Synonyms**

H3/A; H3C2; H3C3; H3C4; H3C6; H3C7; H3C8; H3FA; H3C10; H3C11; H3C12; HIST1H3A; Acetyl-Histone H3-K4/K9/K14/K18/K23/K27


## **Product Information**

Source Isotype **Purification** Rabbit IgG Affinity purification


#### Storage

Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS with 0.05% proclin300,0.05% BSA,50% glycerol,pH7.3.

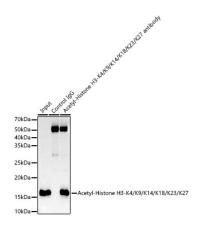

Immunoprecipitation analysis of 600  $\mu$ g extracts of NIH/3T3 cells using 5  $\mu$ g Acetyl-Histone H3-K4/K9/K14/K18/K23/K27 antibody (A22054). Western blot was performed from the immunoprecipitate using Acetyl-Histone H3-K4/K9/K14/K18/K23/K27 antibody (A22054) at a dilution of 1:10000.

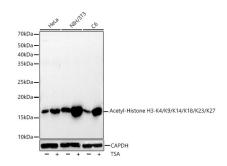


Western blot analysis of various lysates, using Acetyl-Histone H3-K4/K9/K14/K18/K23/K27 Rabbit mAb (A22054) at1:20000 dilution. HeLa,NIH/3T3 and C6 cells were treated by TSA (1 uM) at 37°C for 18 hours. Secondary antibody: HRP-conjugated Goat anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution. Lysates/proteins:  $25\mu g$  per lane.


Blocking buffer: 3% nonfat dry milk in TBST.

Detection: ECL Basic Kit (RM00020).


Exposure time: 30s.

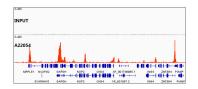


Chromatin immunoprecipitations were performed with cross-linked chromatin from HeLa cells and Acetyl-Histone H3-K4/K9/K14/K18/K23/K27 Rabbit mAb (A22054). The ChIP sequencing results indicate the enrichment pattern of Acetyl-Histone H3-K4/K9/K14/K18/K23/K27 in selected genomic region and representative gene loci (GAPDH), as shown in figure.

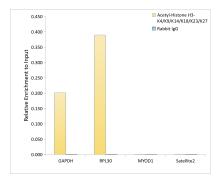


Chromatin immunoprecipitation analysis of extracts of HeLa cells, using Acetyl-Histone H3-K4/K9/K14/K18/K23/K27 antibody (A22054) and rabbit IgG.The amount of immunoprecipitated DNA was checked by quantitative PCR. Histogram was constructed by the ratios of the immunoprecipitated DNA to the input.

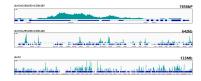


CUT&Tag was performed using the CUT&Tag Assay Kit (pAG-Tn5) for Illumina(RK20265) from  $10^5$  K562 cells with 1  $\mu$ g Acetyl-Histone H3-K4/K9/K14/K18/K23/K27 antibody (A22054), along with a Goat Anti-Rabbit IgG(H+L). The CUT&Tag results indicate the enrichment pattern of Acetyl-Histone H3-K4/K9/K14/K18/K23/K27 in representative gene loci (RPL30), as shown in figure.

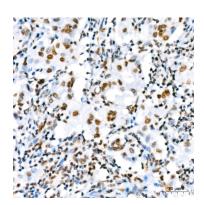


Immunohistochemistry analysis of paraffinembedded Human breast cancer using Acetyl-Histone H3-K4/K9/K14/K18/K23/K27 Rabbit mAb (A22054) at dilution of 1:200 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate Bufferr (pH 6.0) prior to IHC staining.

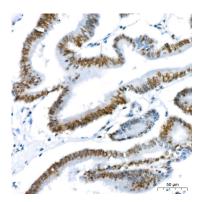


Immunohistochemistry analysis of paraffinembedded Human colon carcinoma using Acetyl-Histone H3-K4/K9/K14/K18/K23/K27 Rabbit mAb (A22054) at dilution of 1:200 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate Bufferr (pH 6.0) prior to IHC staining.

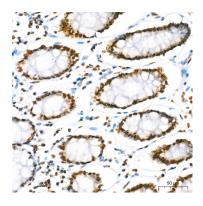


Immunohistochemistry analysis of paraffinembedded Human colon using Acetyl-Histone H3-K4/K9/K14/K18/K23/K27 Rabbit mAb (A22054) at dilution of 1:200 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate Bufferr (pH 6.0) prior to IHC staining.

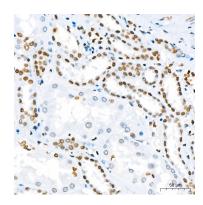


Immunohistochemistry analysis of paraffinembedded Human kidney using Acetyl-Histone H3-K4/K9/K14/K18/K23/K27 Rabbit mAb (A22054) at dilution of 1:200 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate Bufferr (pH 6.0) prior to IHC staining.

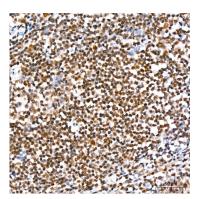


Immunohistochemistry analysis of paraffinembedded Human spleen using Acetyl-Histone H3-K4/K9/K14/K18/K23/K27 Rabbit mAb (A22054) at dilution of 1:200 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate Bufferr (pH 6.0) prior to IHC staining.

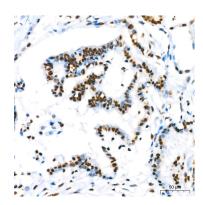

Immunohistochemistry analysis of paraffinembedded Human thyroid cancer using Acetyl-Histone H3-K4/K9/K14/K18/K23/K27 Rabbit mAb (A22054) at dilution of 1:200 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate Bufferr (pH 6.0) prior to IHC staining.

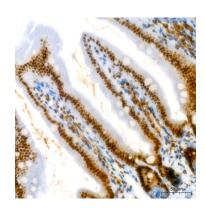

Immunohistochemistry analysis of paraffinembedded Mouse intestin using Acetyl-Histone H3-K4/K9/K14/K18/K23/K27 Rabbit mAb (A22054) at dilution of 1:200 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate Bufferr (pH 6.0) prior to IHC staining.

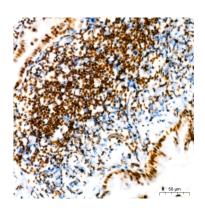

Immunohistochemistry analysis of paraffinembedded Rat lung using Acetyl-Histone H3-K4/K9/K14/K18/K23/K27 Rabbit mAb (A22054) at dilution of 1:200 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate Bufferr (pH 6.0) prior to IHC staining.