# [KO Validated] Non-phospho (Active)β-Catenin -S33/S37/T41 Rabbit mAb

Catalog No.: A22180 KO Validated Recombinant 3 Publications

Basic Information	Background
<b>Observed MW</b> 92kDa	The protein encoded by this gene is part of a complex of proteins that constitute adherens junctions (AJs). AJs are necessary for the creation and maintenance of epithelial cell layers by regulating cell growth and adhesion between cells. The encoded protein also anchors the
Calculated MW 85kDa	actin cytoskeleton and may be responsible for transmitting the contact inhibition signal that causes cells to stop dividing once the epithelial sheet is complete. Finally, this protein binds to the product of the APC gene, which is mutated in adenomatous polyposis of the colon. Mutations in this gene are a cause of colorectal cancer (CRC), pilomatrixoma (PTR), medulloblastoma (MDB), and ovarian cancer. Alternative splicing results in multiple transcript variants.
<b>Category</b> Primary antibody	
Applications WB,IHC-P,ELISA	
<b>Cross-Reactivity</b> Human, Mouse, Rat	
CloneNo number	

## **Recommended Dilutions**

WB	1:2000 - 1:10000
IHC-P	1:1000 - 1:4000
ELISA	Recommended starting concentration is 1 µg/mL. Please optimize the concentration based on your specific assay requirements.

# **Immunogen Information**

Gene ID 1499

**Swiss Prot** P35222

#### Immunogen

Synthetic peptide. This information is considered to be commercially sensitive.

#### **Synonyms**

EVR7; CTNNB; MRD19; NEDSDV; armadillo; Non-phospho (Active)β-Catenin -S33/S37/T41

## Contact

ARC53641

6	400-999-6126
$\mathbf{X}$	cn.market@abclonal.com.cn
€	www.abclonal.com.cn

# **Product Information**

Source Rabbit


Isotype lgG


**Purification** Affinity purification

#### Storage


Store at -20°C. Avoid freeze / thaw cycles. Buffer: PBS with 0.09% Sodium azide, 0.05% BSA, 50% glycerol, pH7.3.



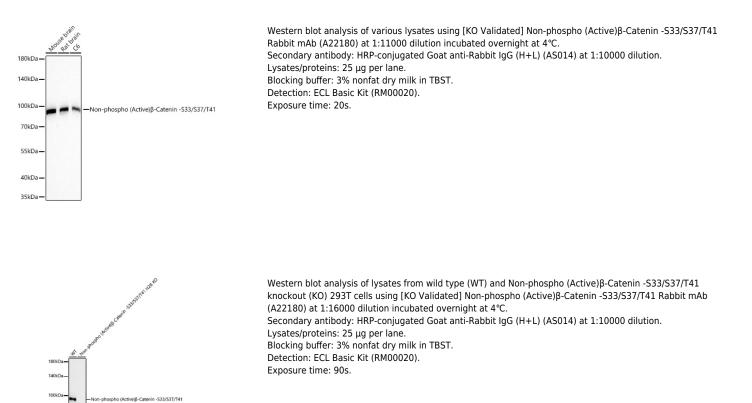


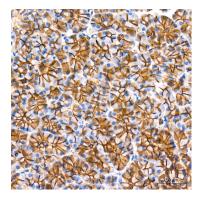

Immunohistochemistry analysis of paraffinembedded Human pancreas tissue using [KO Validated] Non-phospho (Active) $\beta$ -Catenin - S33/S37/T41 Rabbit mAb (A22180) at a dilution of 1:1500 (40x lens). High pressure antigen retrieval performed with 0.01M TrisEDTA Buffer (pH 9.0) prior to IHC staining.

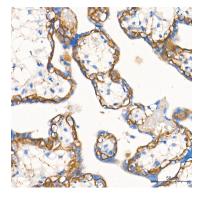


Immunohistochemistry analysis of paraffinembedded Human placenta tissue using [KO Validated] Non-phospho (Active) $\beta$ -Catenin - S33/S37/T41 Rabbit mAb (A22180) at a dilution of 1:1500 (40x lens). High pressure antigen retrieval performed with 0.01M TrisEDTA Buffer (pH 9.0) prior to IHC staining.

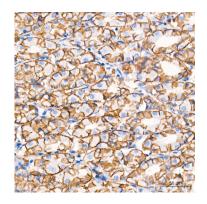


Immunohistochemistry analysis of paraffinembedded Human stomach tissue using [KO Validated] Non-phospho (Active) $\beta$ -Catenin - S33/S37/T41 Rabbit mAb (A22180) at a dilution of 1:1500 (40x lens). High pressure antigen retrieval performed with 0.01M TrisEDTA Buffer (pH 9.0) prior to IHC staining.


### Validation Data

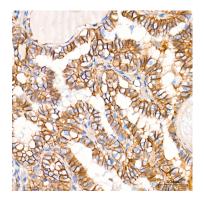


Immunohistochemistry analysis of paraffinembedded Human thyroid cancer tissue using [KO Validated] Non-phospho (Active) $\beta$ -Catenin -S33/S37/T41 Rabbit mAb (A22180) at a dilution of 1:1500 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer (pH 9.0) prior to IHC staining.

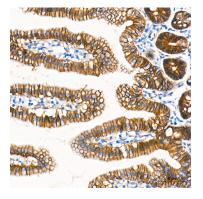

Immunohistochemistry analysis of paraffinembedded Mouse intestin tissue using [KO Validated] Non-phospho (Active) $\beta$ -Catenin - S33/S37/T41 Rabbit mAb (A22180) at a dilution of 1:1500 (40x lens). High pressure antigen retrieval performed with 0.01M TrisEDTA Buffer (pH 9.0) prior to IHC staining.

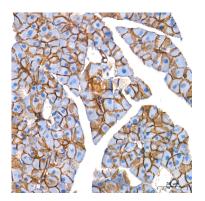

Immunohistochemistry analysis of paraffinembedded Mouse pancreas tissue using [KO Validated] Non-phospho (Active) $\beta$ -Catenin - S33/S37/T41 Rabbit mAb (A22180) at a dilution of 1:1500 (40x lens). High pressure antigen retrieval performed with 0.01M TrisEDTA Buffer (pH 9.0) prior to IHC staining.

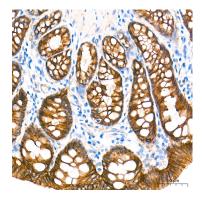

Immunohistochemistry analysis of paraffinembedded Rat colon tissue using [KO Validated] Non-phospho (Active) $\beta$ -Catenin - S33/S37/T41 Rabbit mAb (A22180) at a dilution of 1:1500 (40x lens). High pressure antigen retrieval performed with 0.01M TrisEDTA Buffer (pH 9.0) prior to IHC staining.