Leader in Biomolecular Solutions for Life Science

AHR Rabbit PolymAb®

Catalog No.: A22464PM



Basic Information

Observed MW 105kDa

Calculated MW 96kDa

Category Primary antibody

Applications WB,IHC-P,IF/ICC,ELISA

Cross-Reactivity Human, Mouse

Background

The protein encoded by this gene is a ligand-activated helix-loop-helix transcription factor involved in the regulation of biological responses to planar aromatic hydrocarbons. This receptor has been shown to regulate xenobiotic-metabolizing enzymes such as cytochrome P450. Before ligand binding, the encoded protein is sequestered in the cytoplasm; upon ligand binding, this protein moves to the nucleus and stimulates transcription of target genes.

Recommended Dilutions

WB	1:2000 - 1:10000
IHC-P	1:100 - 1:500
IF/ICC	1:50 - 1:200
ELISA	Recommended starting concentration is 1 µg/mL. Please optimize the concentration based on your specific assay requirements.

Immunogen Information

Gene ID 196

Swiss Prot P35869

Immunogen

Recombinant fusion protein containing a sequence corresponding to amino acids 619-848 of human AHR Rabbit PolymAb $\mbox{@}$ (NP_001612.1).

Synonyms RP85; bHLHe76

Contact


6	400-999-6126
\times	cn.market@abclonal.com.cn
€	www.abclonal.com.cn

Product Information

Source Rabbit **lsotype** IgG Purification Affinity purification

Storage


Store at -20°C. Avoid freeze / thaw cycles. Buffer: PBS with 0.05% proclin300,0.05% BSA,50% glycerol,pH7.3.

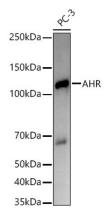


Western blot analysis of lysates from PC-3 cells, using AHR Rabbit PolymAb $\ensuremath{\circledast}$ (A22464PM) at 1:10000 dilution.

Secondary antibody: HRP-conjugated Goat anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution. Lysates/proteins: 25µg per lane. Blocking buffer: 3% nonfat dry milk in TBST. Detection: ECL Basic Kit (RM00020).

Exposure time: 180s.

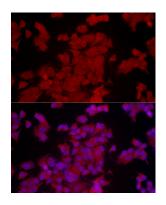


Immunofluorescence analysis of MCF7 cells using AHR Rabbit PolymAb® (A22464PM) at a dilution of 1:100 (40x lens). Secondary antibody: Cy3 Goat anti-Rabbit IgG (H+L) (AS007) at 1:500 dilution. Blue: DAPI for nuclear staining.

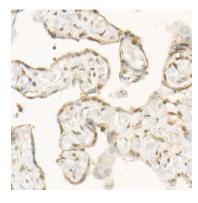


Immunohistochemistry analysis of paraffinembedded Human placenta tissue using AHR Rabbit PolymAb® (A22464PM) at a dilution of 1:500 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate Buffer (pH 6.0) prior to IHC staining.




Immunohistochemistry analysis of paraffinembedded Mouse spleen tissue using AHR Rabbit PolymAb® (A22464PM) at a dilution of 1:500 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate Buffer (pH 6.0) prior to IHC staining.

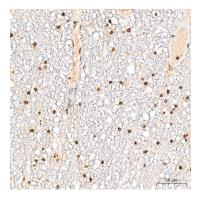


Immunohistochemistry analysis of paraffinembedded Human brain tissue using AHR Rabbit PolymAb® (A22464PM) at a dilution of 1:200 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate Buffer (pH 6.0) prior to IHC staining.

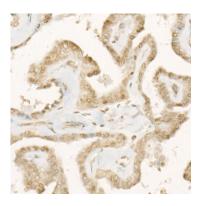


Immunohistochemistry analysis of paraffinembedded Human thyroid cancer tissue using AHR Rabbit PolymAb® (A22464PM) at a dilution of 1:500 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate Buffer (pH 6.0) prior to IHC staining.

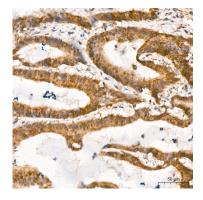


Immunohistochemistry analysis of paraffinembedded Human colon carcinoma tissue using AHR Rabbit PolymAb® (A22464PM) at a dilution of 1:200 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate Buffer (pH 6.0) prior to IHC staining.

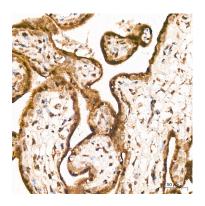


Immunohistochemistry analysis of paraffinembedded Human placenta tissue using AHR Rabbit PolymAb® (A22464PM) at a dilution of 1:200 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate Buffer (pH 6.0) prior to IHC staining.

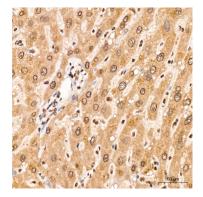


Immunohistochemistry analysis of paraffinembedded Human liver tissue using AHR Rabbit PolymAb® (A22464PM) at a dilution of 1:200 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate Buffer (pH 6.0) prior to IHC staining.

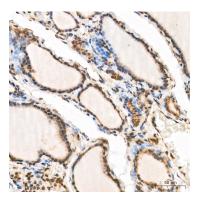

Immunohistochemistry analysis of paraffinembedded Human thyroid tissue using AHR Rabbit PolymAb® (A22464PM) at a dilution of 1:200 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate Buffer (pH 6.0) prior to IHC staining.

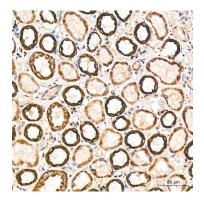

Immunohistochemistry analysis of paraffinembedded Human kidney tissue using AHR Rabbit PolymAb® (A22464PM) at a dilution of 1:200 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate Buffer (pH 6.0) prior to IHC staining.

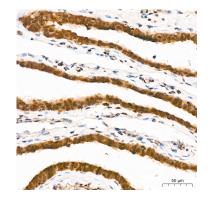

Immunohistochemistry analysis of paraffinembedded Human fallopian tube tissue using AHR Rabbit PolymAb® (A22464PM) at a dilution of 1:200 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate Buffer (pH 6.0) prior to IHC staining.