

ABflo® 488 Rabbit anti-Human Ig κ Light Chain mAb

Catalog No.: A23900

Basic Information

Observed MW

Calculated MW
11KDa

Category
Primary antibody

Applications
FC

Cross-Reactivity
Human

CloneNo number
ARC60646-ABflo488

Conjugate
ABflo® 488. Ex:491nm. Em:516nm.

Recommended Dilutions

FC 5 μl per 10⁶ cells in
100 μl volume

Contact

☎ | 400-999-6126

✉ | cn.market@abclonal.com.cn

🌐 | www.abclonal.com.cn

Background

Immunoglobulins belong to a group of related glyco proteins which make up 20% of serum proteins. Antigens and immunoglobulins react to confer immunity to individuals. Immunoglobulins have similar structures of two identical heavy chains and two identical light chains. Both the heavy chains and the light chains are divided into constant and variable regions. The constant regions have the same amino acid sequences between all the immunoglobulin classes. The variable regions have approximately 110 amino acids with high sequence variability. The amino acid sequence of the heavy chain determines the class of an immunoglobulin. The five types of immunoglobulin heavy chains are known as: IgG, IgA, IgM, IgD, and IgE. IgG is divided into four subclasses, and IgA is divided into two subclasses. In serum IgA and IgG are monomers with a single 4 polypeptide unit; while, IgM is a pentamer. IgA may also form polymers. Kappa light chain antibody can be used for the identification of leukemias, plasmacytomas and certain non Hodgkin's Lymphomas. Kappa light chain contains one immunoglobulin like domain. The EU sequence has the INV allotypic marker, Ala 45 and Val 83. The ROY sequence has the INV allotypic marker, Ala 45 and Leu 83.

Immunogen Information

Gene ID
3514

Swiss Prot
P01834

Immunogen

Recombinant fusion protein containing a sequence corresponding to amino acids 1-107 of human Ig κ Light Chain(P01834).

Synonyms

IGKC;HCAK1; IGKCD; Km

Product Information

Source
Rabbit

Isotype
IgG


Purification
Affinity purification

Storage

Store at 2-8°C. Avoid freeze.

Buffer: PBS with 0.03% proclin300,0.2% BSA,pH7.3.

Validation Data

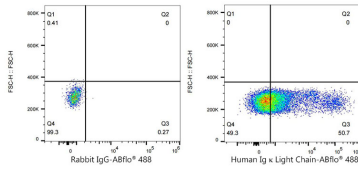

Flow cytometry: 1×10^6 Human PBMC were surface-stained with ABflo® 488 Rabbit IgG isotype control (A22069,5 μ l/Test, left) or ABflo® 488 Rabbit anti-Human Ig κ Light Chain mAb (ARC60646-01-01-ABf488,5 μ l/Test, right).