Leader in Biomolecular Solutions for Life Science

ROCK2 Rabbit mAb

Catalog No.: A2395 Recombinant 8 Publications



Basic Information

Observed MW 161kDa

Calculated MW 161kDa

Category Primary antibody

Applications WB,IHC-P,ELISA

Cross-Reactivity Human, Mouse, Rat

CloneNo number ARC0744

Background

The protein encoded by this gene is a serine/threonine kinase that regulates cytokinesis, smooth muscle contraction, the formation of actin stress fibers and focal adhesions, and the activation of the c-fos serum response element. This protein, which is an isozyme of ROCK1 is a target for the small GTPase Rho.

Recommended Dilutions

WB	1:1000 - 1:2000
IHC-P	1:200 - 1:800
ELISA	Recommended starting concentration is 1 µg/mL. Please optimize the concentration based on your specific assay requirements.

Immunogen Information

Gene ID 9475

Swiss Prot 075116

Immunogen

Synthetic peptide. This information is considered to be commercially sensitive.

Synonyms ROCK-II; ROCK2

Contact

6	400-999-6126
\times	cn.market@abclonal.com.cn
€	www.abclonal.com.cn

Product Information

Source Rabbit


Isotype lgG

Purification Affinity purification


Storage

Store at -20°C. Avoid freeze / thaw cycles. Buffer: PBS with 0.02% sodium azide, 0.05% BSA, 50% glycerol, pH7.3.


Validation Data

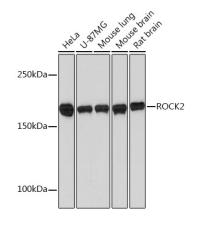

Western blot analysis of various lysates using ROCK2 Rabbit mAb (A2395) at 1[]1000 dilution. Secondary antibody: HRP-conjugated Goat anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution. Lysates/proteins: 25µg per lane. Blocking buffer: 3% nonfat dry milk in TBST. Detection: ECL Basic Kit (RM00020). Exposure time: 90s.

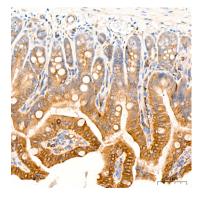

Immunohistochemistry analysis of paraffinembedded Mouse intestin tissue using ROCK2 Rabbit mAb (A2395) at a dilution of 1:200 (40x lens). High pressure antigen retrieval was performed with 0.01 M citrate buffer (pH 6.0) prior to IHC staining.

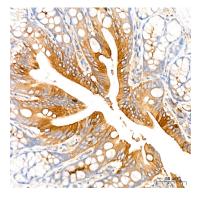

Immunohistochemistry analysis of paraffinembedded Rat colon tissue using ROCK2 Rabbit mAb (A2395) at a dilution of 1:200 (40x lens). High pressure antigen retrieval was performed with 0.01 M citrate buffer (pH 6.0) prior to IHC staining.