# Collagen I/COL1A1 Rabbit PolymAb®

www.abclonal.com

ABclonal

Catalog No.: A24112PM

## **Basic Information**

#### **Observed MW**

220kDa

#### **Calculated MW**

139kDa

#### Category

Primary antibody

## **Applications**

WB,IF/ICC,IHC-P,ELISA

#### **Cross-Reactivity**

Human, Mouse, Rat

# **Background**

This gene encodes the pro-alpha1 chains of type I collagen whose triple helix comprises two alpha1 chains and one alpha2 chain. Type I is a fibril-forming collagen found in most connective tissues and is abundant in bone, cornea, dermis and tendon. Mutations in this gene are associated with osteogenesis imperfecta types I-IV, Ehlers-Danlos syndrome type VIIA, Ehlers-Danlos syndrome Classical type, Caffey Disease and idiopathic osteoporosis. Reciprocal translocations between chromosomes 17 and 22, where this gene and the gene for platelet-derived growth factor beta are located, are associated with a particular type of skin tumor called dermatofibrosarcoma protuberans, resulting from unregulated expression of the growth factor. Two transcripts, resulting from the use of alternate polyadenylation signals, have been identified for this gene.

## **Recommended Dilutions**

**WB** 1:6000 - 1:36000

**IF/ICC** 1:200 - 1:2000

IHC-P 1:2000 - 1:8000

**ELISA** Recommended starting

concentration is 1 µg/mL.

Please optimize the
concentration based on
your specific assay
requirements.

# Immunogen Information

**Gene ID**1277

Swiss Prot
P02452

#### **Immunogen**

Synthetic peptide. This information is considered to be commercially sensitive.

## **Synonyms**

OI1; OI2; OI3; OI4; EDSC; CAFYD; EDSARTH1

## **Contact**

2		400-999-6126
$\bowtie$		cn.market@abclonal.com.cn
•	ī	www.abclonal.com.cn


## **Product Information**

SourceIsotypePurificationRabbitIgGAffinity purification

#### Storage

Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS with 0.09% sodium azide, 0.05% BSA, 50% glycerol, pH7.3.

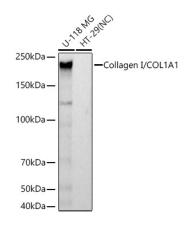

Western blot analysis of various lysates using Collagen I/COL1A1 Rabbit PolymAb® (A24112PM) at 1:6000 dilution incubated overnight at  $4^{\circ}$ C.


Secondary antibody: HRP-conjugated Goat anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.

Lysates/proteins: 25 µg per lane.


Blocking buffer: 3% nonfat dry milk in TBST.

Detection: ECL Basic Kit (RM00020). Negative control (NC): HT-29


Exposure time: 10s.

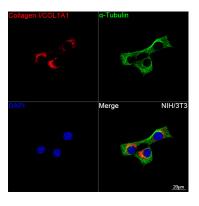


Confocal imaging of NIH/3T3 cells using Collagen I/COL1A1 Rabbit PolymAb (A24112PM, dilution 1:200) followed by a further incubation with Cy3 Goat Anti-Rabbit IgG (H+L) (AS007, dilution 1:500) (Red). The cells were counterstained with  $\alpha$ -Tubulin Mouse mAb (AC012, dilution 1:400) followed by incubation with ABflo (8 488-conjugated Goat Anti-Mouse IgG (H+L) Ab (AS076, dilution 1:500) (Green). DAPI was used for nuclear staining (Blue). Objective: 100x.

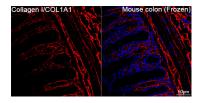


Confocal imaging of frozen sections Mouse colon tissue using Collagen I/COL1A1 Rabbit PolymAb® (A24112PM, dilution 1:200) followed by a further incubation with Cy3 Goat Anti-Rabbit IgG (H+L) (AS007, dilution 1:500) (Red). DAPI was used for nuclear staining (Blue). Microwave antigen retrieval performed with 0.01M Citrate Buffer (pH 6.0) prior to IF staining. Objective: 40x.

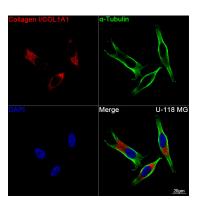


Confocal imaging of U-118 MG cells using Collagen I/COL1A1 Rabbit PolymAb (A24112PM, dilution 1:200) followed by a further incubation with Cy3 Goat Anti-Rabbit IgG (H+L) (AS007, dilution 1:500) (Red). The cells were counterstained with  $\alpha$ -Tubulin Mouse mAb (AC012, dilution 1:400) followed by incubation with ABflo (B48-conjugated Goat Anti-Mouse IgG (H+L) Ab (AS076, dilution 1:500) (Green). DAPI was used for nuclear staining (Blue). Objective: 100x.

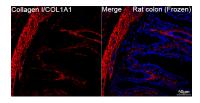


Confocal imaging of frozen sections Rat colon tissue using Collagen I/COL1A1 Rabbit PolymAb® (A24112PM, dilution 1:200) followed by a further incubation with Cy3 Goat Anti-Rabbit IgG (H+L) (AS007, dilution 1:500) (Red). DAPI was used for nuclear staining (Blue). Microwave antigen retrieval performed with 0.01M Citrate Buffer (pH 6.0) prior to IF staining. Objective: 40x.

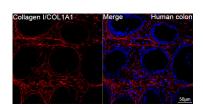


Confocal imaging of paraffin-embedded Human colon tissue using Collagen I/COL1A1 Rabbit PolymAb® (A24112PM, dilution 1:200) followed by a further incubation with Cy3 Goat Anti-Rabbit IgG (H+L) (A5007, dilution 1:500) (Red). DAPI was used for nuclear staining (Blue). High pressure antigen retrieval performed with 0.01M Citrate Buffer (pH 6.0) prior to IF staining. Objective: 40x.

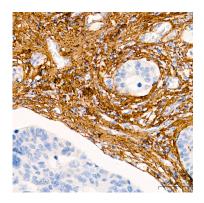


Immunohistochemistry analysis of paraffinembedded Human colon carcinoma tissue using Collagen I/COL1A1 Rabbit PolymAb® (A24112PM) at a dilution of 1:6000 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer (pH 9.0) prior to IHC staining.

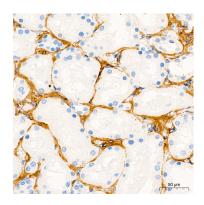

Immunohistochemistry analysis of paraffinembedded Human kidney tissue using Collagen I/COL1A1 Rabbit PolymAb® (A24112PM) at a dilution of 1:6000 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer (pH 9.0) prior to IHC staining.

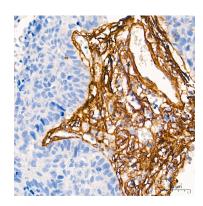


Immunohistochemistry analysis of paraffinembedded Human lung cancer tissue using Collagen I/COL1A1 Rabbit PolymAb® (A24112PM) at a dilution of 1:6000 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer (pH 9.0) prior to IHC staining.

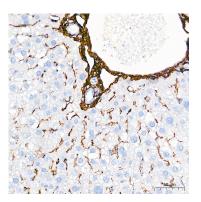

Immunohistochemistry analysis of paraffinembedded Mouse liver tissue using Collagen I/COL1A1 Rabbit PolymAb® (A24112PM) at a dilution of 1:6000 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer (pH 9.0) prior to IHC staining.



Immunohistochemistry analysis of paraffinembedded Mouse lung tissue using Collagen I/COL1A1 Rabbit PolymAb® (A24112PM) at a dilution of 1:6000 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer (pH 9.0) prior to IHC staining.

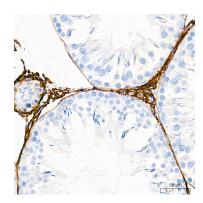

Immunohistochemistry analysis of paraffinembedded Rat testis tissue using Collagen I/COL1A1 Rabbit PolymAb® (A24112PM) at a dilution of 1:6000 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer (pH 9.0) prior to IHC staining.