# **SOX9 Rabbit pAb**

Catalog No.: A2479 8 Publications



## **Basic Information**

#### **Observed MW**

70kDa

### **Calculated MW**

56kDa

### Category

Primary antibody

### **Applications**

ELISA,WB,IHC-P

### **Cross-Reactivity**

Human, Mouse, Rat

## **Background**

The protein encoded by this gene recognizes the sequence CCTTGAG along with other members of the HMG-box class DNA-binding proteins. It acts during chondrocyte differentiation and, with steroidogenic factor 1, regulates transcription of the anti-Muellerian hormone (AMH) gene. Deficiencies lead to the skeletal malformation syndrome campomelic dysplasia, frequently with sex reversal.

## **Recommended Dilutions**

**WB** 1:500 - 1:2000

**IHC-P** 1:50 - 1:200

## **Immunogen Information**

**Gene ID Swiss Prot**6662
P48436

### **Immunogen**

Recombinant fusion protein containing a sequence corresponding to amino acids 250-339 of human SOX9 (NP\_000337.1).

### **Synonyms**

CMD1; SRA1; CMPD1; SRXX2; SRXY10; SOX9

### **Contact**

<u>a</u>		400-999-6126
$\bowtie$		cn.market@abclonal.com.cn
$\odot$	П	www.abclonal.com.cn


### **Product Information**

SourceIsotypePurificationRabbitIgGAffinity purification

### **Storage**

Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS with 0.05% proclin300,50% glycerol,pH7.3.

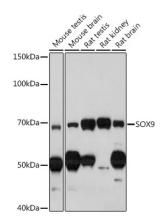


Western blot analysis of various lysates using SOX9 Rabbit pAb (A2479) at 1:1000 dilution. Secondary antibody: HRP Goat Anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.


Lysates/proteins: 25µg per lane.

Blocking buffer: 3% nonfat dry milk in TBST. Detection: ECL Basic Kit (RM00020).


Exposure time: 120s.

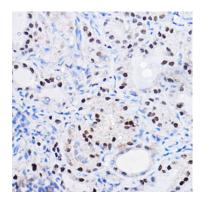

Immunohistochemistry analysis of paraffinembedded human thyroid cancer using SOX9 Rabbit pAb (A2479) at dilution of 1:100 (40x lens).Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.

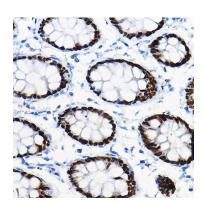

Immunohistochemistry analysis of paraffinembedded human colon using SOX9 Rabbit pAb (A2479) at dilution of 1:100 (40x lens). Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.

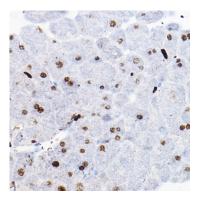

Immunohistochemistry analysis of paraffinembedded mouse pancreas using SOX9 Rabbit pAb (A2479) at dilution of 1:100 (40x lens). Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.