



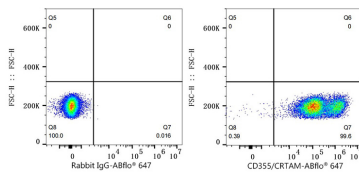

## Validation Data

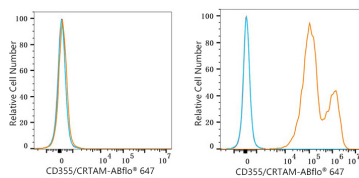

Flow cytometry:  $1 \times 10^6$  293T cells (negative control, left) and 293T (Transfection, right) cells were surface-stained with ABflo® 647 Rabbit anti-Human CD355/CRTAM mAb (A24819,5  $\mu$ l/Test, orange line) or ABflo® 647 Rabbit IgG isotype control (A22070,5  $\mu$ l/Test, blue line). Non-fluorescently stained cells were used as blank control (red line).

Flow cytometry:  $1 \times 10^6$  293T (Transfection) cells were surface-stained with ABflo® 647 Rabbit IgG isotype control (A22070,5  $\mu$ l/Test, left) or ABflo® 647 Rabbit anti-Human CD355/CRTAM mAb (A24819,5  $\mu$ l/Test, right).