


Validation Data

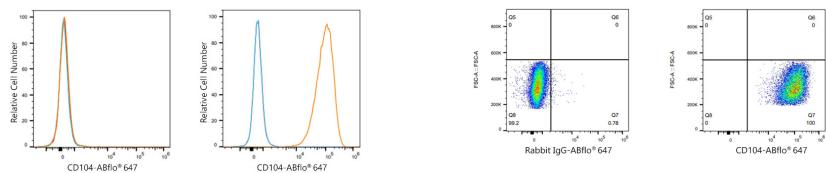

Flow cytometry: 1X10⁶ K-562 cells (negative control, left) and A431 cells (right) were surface-stained with CD104/Integrin β 4 Rabbit mAb (A25086, 2 μ g/mL, orange line) or ABflo® 647 Rabbit IgG isotype control (A22070, 5 μ l/Test, blue line), followed by Alexa Fluor® 647 conjugated goat anti-rabbit pAb staining. Non-fluorescently stained cells were used as blank control (red line).

Flow cytometry: 1X10⁶ A431 cells were surface-stained with ABflo® 647 Rabbit IgG isotype control (A22070, 5 μ l/Test, left) or CD104/Integrin β 4 Rabbit mAb (A25086, 2 μ g/mL, right).