# PIK3CG Rabbit pAb

Catalog No.: A25468



# **Basic Information**

#### **Observed MW**

Refer to figures

### **Calculated MW**

126kDa

### Category

Primary antibody

### **Applications**

IHC-P,ELISA

#### **Cross-Reactivity**

Human

# **Background**

Phosphoinositide 3-kinases (PI3Ks) phosphorylate inositol lipids and are involved in the immune response. The protein encoded by this gene is a class I catalytic subunit of PI3K. Like other class I catalytic subunits (p110-alpha p110-beta, and p110-delta), the encoded protein binds a p85 regulatory subunit to form PI3K. This gene is located in a commonly deleted segment of chromosome 7 previously identified in myeloid leukemias. Several transcript variants encoding the same protein have been found for this gene.

# **Recommended Dilutions**

**IHC-P** 

1:50 - 1:200

**ELISA** 

Recommended starting concentration is 1 µg/mL.
Please optimize the concentration based on your specific assay requirements.

# Immunogen Information

Gene ID

Swiss Prot P48736

5294

#### **Immunogen**

Recombinant protein (or fragment). This information is considered to be commercially sensitive.

# **Synonyms**

PI3K; PIK3; IMD97; PI3CG; PI3Kgamma; p110gamma; p120-PI3K; PIK3CG

# **Contact**

	400-999-6126
<b>×</b>	cn.market@abclonal.com.cn
	www.abclonal.com.cn

# **Product Information**


SourceIsotypePurificationRabbitIgGAffinity purification

#### Storage


Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS with 0.05% proclin300,50% glycerol,pH7.3.

# **Validation Data**

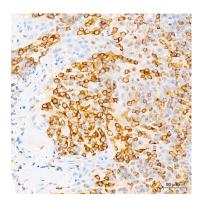

Immunohistochemistry analysis of paraffinembedded Human tonsil tissue using PI 3 Kinase catalytic subunit gamma Rabbit pAb (A25468) at a dilution of 1:1200 (40x lens). High pressure antigen retrieval was performed with 0.01 M citrate buffer (pH 6.0) prior to IHC staining.

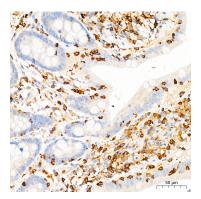

Immunohistochemistry analysis of paraffinembedded Human colon tissue using PI 3 Kinase catalytic subunit gamma Rabbit pAb (A25468) at a dilution of 1:1200 (40x lens). High pressure antigen retrieval was performed with 0.01 M citrate buffer (pH 6.0) prior to IHC staining.