Insulin Rabbit pAb

Catalog No.: A25589



Basic Information

Observed MW

Refer to figures

Calculated MW

12kDa

Category

Primary antibody

Applications

IF-P,IHC-P,ELISA

Cross-Reactivity

Human

Background

This gene encodes insulin, a peptide hormone that plays a vital role in the regulation of carbohydrate and lipid metabolism. After removal of the precursor signal peptide, proinsulin is post-translationally cleaved into three peptides: the B chain and A chain peptides, which are covalently linked via two disulfide bonds to form insulin, and C-peptide. Binding of insulin to the insulin receptor (INSR) stimulates glucose uptake. A multitude of mutant alleles with phenotypic effects have been identified, including insulin-dependent diabetes mellitus, permanent neonatal diabetes diabetes mellitus, maturity-onset diabetes of the young type 10 and hyperproinsulinemia. There is a read-through gene, INS-IGF2, which overlaps with this gene at the 5' region and with the IGF2 gene at the 3' region.

Recommended Dilutions

IF-P 1:50 - 1:200

IHC-P 1:50 - 1:200

 $\begin{array}{cc} \textbf{ELISA} & \text{Recommended starting} \\ & \text{concentration is 1 } \mu\text{g/mL}. \end{array}$

Please optimize the concentration based on your specific assay requirements.

Immunogen Information

Gene ID3630

Swiss Prot
P01308

Immunogen

Synthetic peptide. This information is considered to be commercially sensitive.

Synonyms

IDDM; ILPR; IRDN; IDDM1; IDDM2; PNDM4; MODY10

Contact

a		400-999-6126
\bowtie		cn.market@abclonal.com.cn
•	Ī	www.abclonal.com.cn

Product Information


SourceIsotypePurificationRabbitIgGAffinity purification


Storage


Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS with 0.09% Sodium azide,50% glycerol,pH7.3.

Validation Data

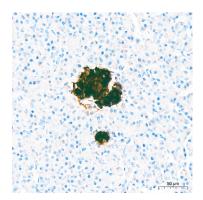

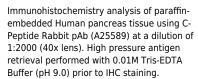

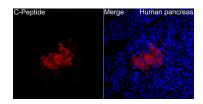

Immunofluorescence analysis of Human pancreas tissue using C-Peptide Rabbit pAb (A25589) at a dilution of 1:100 (40x lens). Secondary antibody: Cy3-conjugated Goat anti-Rabbit IgG (H+L)(AS007) at 1:500 dilution. Blue: DAPI for nuclear staining. High pressure antigen retrieval performed with 0.01M Citrate Buffer (pH 6.0) prior to IF staining.