


Validation Data

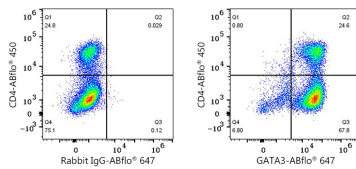

Flow cytometry: 1×10^6 Human PBMC (treated with 50ng/ml PMA + 1 μ g/ml Ionomycin + 2 μ M Monensin for 6 hours) were surface-stained with ABflo® 450 Rabbit anti-Human/Monkey CD4 mAb (A26523, 5 μ l/Test), then intracellularly-stained ABflo® 647 Rabbit IgG isotype control (A22070, 5 μ l/Test), left) or ABflo® 647 Rabbit anti-Human GATA3 mAb (A25946, 5 μ l/Test, right). Cells in the lymphocyte gate were used for analysis.