

[KD Validated] VDAC3 Rabbit mAb

Catalog No.: A26244 **Recombinant**

Basic Information

Observed MW

31kDa

Calculated MW

31kDa

Category

Primary antibody

Applications

WB,IHC-P,IF/ICC,ELISA

Cross-Reactivity

Human, Mouse, Rat

CloneNo number

ARC69554

Background

This gene encodes a voltage-dependent anion channel (VDAC), and belongs to the mitochondrial porin family. VDACS are small, integral membrane proteins that traverse the outer mitochondrial membrane and conduct ATP and other small metabolites. They are known to bind several kinases of intermediary metabolism, thought to be involved in translocation of adenine nucleotides, and are hypothesized to form part of the mitochondrial permeability transition pore, which results in the release of cytochrome c at the onset of apoptotic cell death. Alternatively transcript variants encoding different isoforms have been described for this gene.

Recommended Dilutions

WB	1:500 - 1:1000
IHC-P	1:50 - 1:200
IF/ICC	1:50 - 1:200
ELISA	Recommended starting concentration is 1 µg/mL. Please optimize the concentration based on your specific assay requirements.

Contact

	400-999-6126
	cn.market@abclonal.com.cn
	www.abclonal.com.cn

Immunogen Information

Gene ID

7419

Swiss Prot

Q9Y277

Immunogen

A synthetic peptide corresponding to a sequence within amino acids 1-100 of human VDAC3 (NP_005653.3).

Synonyms

VDAC-3; HD-VDAC3

Product Information

Source

Rabbit

Isotype

IgG

Purification


Affinity purification

Storage


Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS with 0.09% Sodium azide,0.05% BSA,50% glycerol,pH7.3.


Validation Data

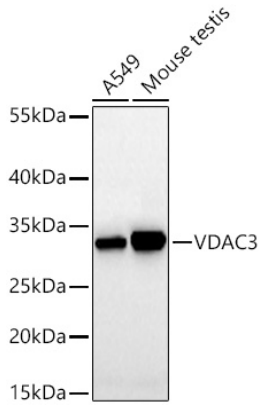

Western blot analysis of various lysates using [KD Validated] VDAC3 Rabbit mAb (A26244) at 1:1000 dilution incubated overnight at 4°C.
Secondary antibody: Goat anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.
Lysates/proteins: 25 µg per lane.
Blocking buffer: 3% nonfat dry milk in TBST.
Detection: ECL Basic Kit (RM00020).
Exposure time: 20s.

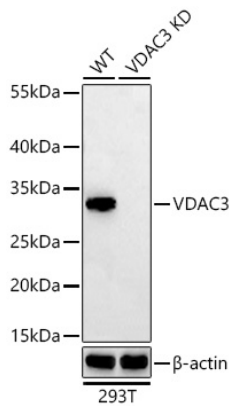


Western blot analysis of lysates from wild type (WT) and VDAC3 knockdown (KD) 293T cells using [KD Validated] VDAC3 Rabbit mAb (A26244) at 1:1000 dilution incubated overnight at 4°C.
Secondary antibody: HRP-conjugated Goat anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.
Lysates/proteins: 25 µg per lane.
Blocking buffer: 3% nonfat dry milk in TBST.
Detection: ECL Basic Kit (RM00020).
Exposure time: 20s.

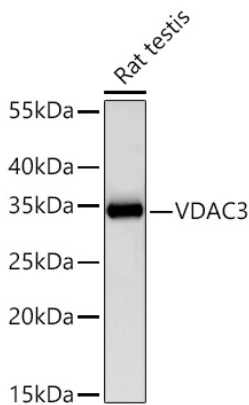


Western blot analysis of lysates from Rat testis using [KD Validated] VDAC3 Rabbit mAb (A26244) at 1:1000 dilution incubated overnight at 4°C.
Secondary antibody: HRP-conjugated Goat anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.
Lysates/proteins: 25 µg per lane.
Blocking buffer: 3% nonfat dry milk in TBST.
Detection: ECL Basic Kit (RM00020).
Exposure time: 15s.


Validation Data

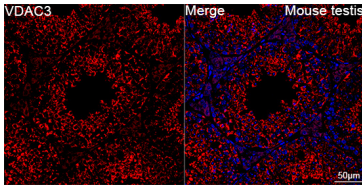


Confocal imaging of paraffin-embedded Mouse testis tissue using [KD Validated] VDAC3 Rabbit mAb (A26244, dilution 1:200) followed by a further incubation with Cy3 Goat Anti-Rabbit IgG (H+L) (AS007, dilution 1:500) (Red). DAPI was used for nuclear staining (Blue). High pressure antigen retrieval performed with 0.01M Citrate Buffer(pH 6.0) prior to IF staining. Objective: 40x.

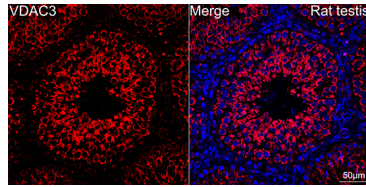


Confocal imaging of paraffin-embedded Rat testis tissue using [KD Validated] VDAC3 Rabbit mAb (A26244, dilution 1:200) followed by a further incubation with Cy3 Goat Anti-Rabbit IgG (H+L) (AS007, dilution 1:500) (Red). DAPI was used for nuclear staining (Blue). High pressure antigen retrieval performed with 0.01M Citrate Buffer(pH 6.0) prior to IF staining. Objective: 40x.

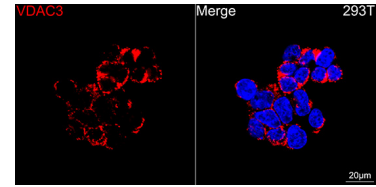

Confocal imaging of 293T cells using [KD Validated] VDAC3 Rabbit mAb (A26244, dilution 1:200) followed by a further incubation with Cy3 Goat Anti-Rabbit IgG (H+L) (AS007, dilution 1:500) (Red). DAPI was used for nuclear staining (Blue). Objective: 100x.

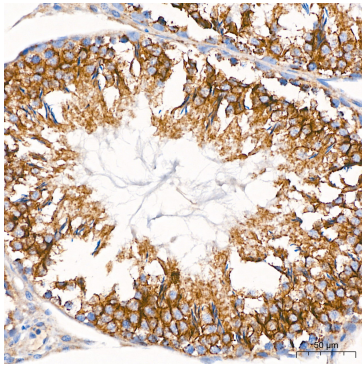

Immunohistochemistry analysis of paraffin-embedded Rat testis tissue using [KD Validated] VDAC3 Rabbit mAb (A26244) at a dilution of 1:200 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate Buffer(pH 6.0) prior to IHC staining.

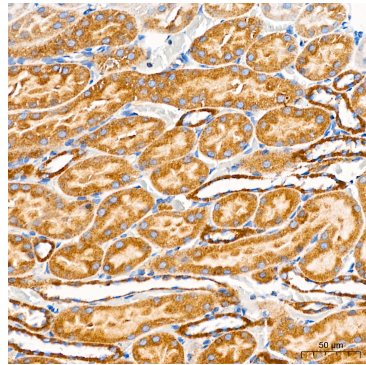

Immunohistochemistry analysis of paraffin-embedded Rat kidney tissue using [KD Validated] VDAC3 Rabbit mAb (A26244) at a dilution of 1:200 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate Buffer(pH 6.0) prior to IHC staining.