
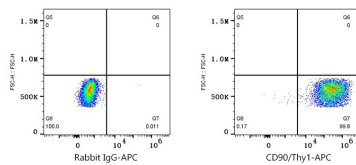


Validation Data

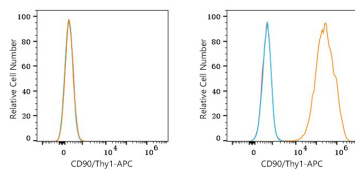

Flow cytometry: 1×10^6 K-562 cells (negative control, left) and U-251 MG cells (right) were surface-stained with APC Rabbit anti-Human CD90/Thy1 mAb (A26399, 5 μ l/Test, orange line) or APC Rabbit IgG isotype control (A24173, 5 μ l/Test, blue line). Non-fluorescently stained cells were used as blank control (red line).

Flow cytometry: 1×10^6 U-251 MG cells were surface-stained with APC Rabbit IgG isotype control (A24173, 5 μ l/Test, left) or APC Rabbit anti-Human CD90/Thy1 mAb (A26399, 5 μ l/Test, right).