# **GOLPH2 Rabbit PolymAb®**

Catalog No.: A26409PM



### **Basic Information**

#### **Observed MW**

80kDa

#### **Calculated MW**

45kDa

### Category

Primary antibody

#### **Applications**

WB,IHC-P,IF/ICC,ELISA

#### **Cross-Reactivity**

Human, Mouse, Rat

## **Background**

The Golgi complex plays a key role in the sorting and modification of proteins exported from the endoplasmic reticulum. The protein encoded by this gene is a type II Golgi transmembrane protein. It processes proteins synthesized in the rough endoplasmic reticulum and assists in the transport of protein cargo through the Golgi apparatus. The expression of this gene has been observed to be upregulated in response to viral infection. Alternatively spliced transcript variants encoding the same protein have been described for this gene.

### **Recommended Dilutions**

**WB** 1:1000 - 1:3000

IHC-P 1:5000 - 1:20000

**IF/ICC** 1:1000 - 1:4000

**ELISA** Recommended starting

concentration is 1 µg/mL.

Please optimize the
concentration based on
your specific assay
requirements.

# **Immunogen Information**

**Gene ID**51280

Swiss Prot
Q8NB|4

#### **Immunogen**

Recombinant fusion protein containing a sequence corresponding to amino acids 250-350 of human GOLPH2. (NP\_057632.2).

### Synonyms

GP73; HEL46; GOLPH2; C9orf155; PSEC0257; bA379P1.3

### **Contact**

<b>a</b>		400-999-6126
$\bowtie$		cn.market@abclonal.com.cn
$\odot$	Ī	www.abclonal.com.cn


### **Product Information**

SourceIsotypePurificationRabbitIgGAffinity purification

#### Storage

Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS with 0.09% Sodium azide, 0.05% BSA, 50% glycerol, pH7.3.

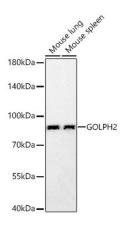

Western blot analysis of various lysates using GOLPH2 Rabbit PolymAb® (A26409PM) at 1:3000 dilution incubated overnight at  $4^{\circ}$ C.


Secondary antibody: HRP-conjugated Goat anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.

Lysates/proteins: 25 µg per lane.


Blocking buffer: 3% nonfat dry milk in TBST.

Detection: ECL Basic Kit (RM00020).


Exposure time: 60s.

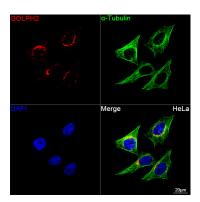


Confocal imaging of HeLa cells using GOLPH2 Rabbit PolymAb® (A26409PM, dilution 1:2000) followed by a further incubation with Cy3 Goat Anti-Rabbit IgG (H+L) (AS007, dilution 1:500) (Red). The cells were counterstained with  $\alpha$ -Tubulin Mouse mAb (AC012, dilution 1:400) followed by incubation with ABflo® 488-conjugated Goat Anti-Mouse IgG (H+L) Ab (AS076, dilution 1:500) (Green). DAPI was used for nuclear staining (Blue). Objective: 100x.

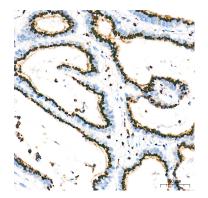


Immunohistochemistry analysis of paraffinembedded Human colon carcinoma tissue using GOLPH2 Rabbit PolymAb® (A26409PM) at a dilution of 1:8000 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer(pH 9.0) prior to IHC staining.

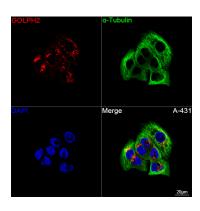


Confocal imaging of A-431 cells using GOLPH2 Rabbit PolymAb® (A26409PM, dilution 1:2000) followed by a further incubation with Cy3 Goat Anti-Rabbit IgG (H+L) (AS007, dilution 1:500) (Red). The cells were counterstained with  $\alpha$ -Tubulin Mouse mAb (AC012, dilution 1:400) followed by incubation with ABflo® 488-conjugated Goat Anti-Mouse IgG (H+L) Ab (AS076, dilution 1:500) (Green). DAPI was used for nuclear staining (Blue). Objective: 100x.

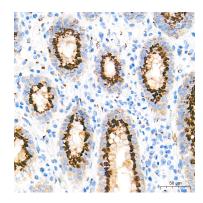


Immunohistochemistry analysis of paraffinembedded Human colon tissue using GOLPH2 Rabbit PolymAb® (A26409PM) at a dilution of 1:8000 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer(pH 9.0) prior to IHC staining.

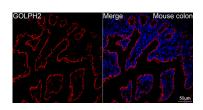


Confocal imaging of paraffin-embedded Mouse colon tissue using GOLPH2 Rabbit PolymAb® (A26409PM, dilution 1:2000) followed by a further incubation with Cy3 Goat Anti-Rabbit IgG (H+L) (AS007, dilution 1:500) (Red). DAPI was used for nuclear staining (Blue). High pressure antigen retrieval performed with 0.01M Citrate Buffer (pH 6.0) prior to IF staining. Objective: 40x.

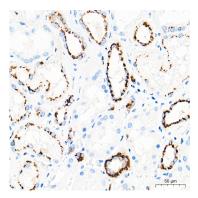

Immunohistochemistry analysis of paraffinembedded Human kidney tissue using GOLPH2 Rabbit PolymAb® (A26409PM) at a dilution of 1:8000 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer(pH 9.0) prior to IHC staining.

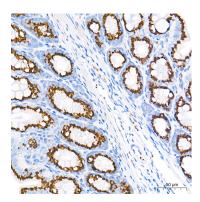

Immunohistochemistry analysis of paraffinembedded Mouse intestin tissue using GOLPH2 Rabbit PolymAb® (A26409PM) at a dilution of 1:8000 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer(pH 9.0) prior to IHC staining.

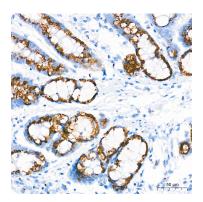

Immunohistochemistry analysis of paraffinembedded Rat colon tissue using GOLPH2 Rabbit PolymAb® (A26409PM) at a dilution of 1:8000 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer(pH 9.0) prior to IHC staining.