ICAM-1/CD54 Rabbit PolymAb®

Catalog No.: A26412-PM



Basic Information

Observed MW

82-100kDa

Calculated MW

58kDa/59kDa

Category

Primary antibody

Applications

WB,IHC-P,IF/ICC,ELISA

Cross-Reactivity

Human, Mouse, Rat

Background

This gene encodes a cell surface glycoprotein which is typically expressed on endothelial cells and cells of the immune system. It binds to integrins of type CD11a / CD18, or CD11b / CD18 and is also exploited by Rhinovirus as a receptor.

Recommended Dilutions

WB 1:1000 - 1:30000

IHC-P 1:2000 - 1:10000

IF/ICC 1:2000 - 1:4000

ELISA Recommended starting

concentration is 1 µg/mL.

Please optimize the concentration based on your specific assay requirements.

Immunogen Information

 Gene ID
 Swiss Prot

 3383/15894
 P05362/P13597

Immunogen

Recombinant protein of human/mouse ICAM-1/CD54.

Synonyms

MAL; CD54; Icam; Ly-4; Ly-47; Icam-1; MALA-2

Contact

a	400-999-6126
×	cn.market@abclonal.com.cn
$\overline{\Box}$	www.ahclonal.com.cn


Product Information

SourceIsotypePurificationRabbitIgGAffinity purification

Storage

Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS with 0.09% Sodium azide, 0.05% BSA, 50% glycerol, pH7.3.

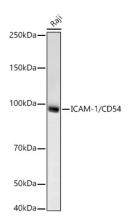

Western blot analysis of lysates from Raji cells using ICAM-1/CD54 Rabbit PolymAb® (A26412-PM) at 1:1000 dilution incubated overnight at 4°C.


Secondary antibody: HRP-conjugated Goat anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.

Lysates/proteins: 25 µg per lane.

Blocking buffer: 3% nonfat dry milk in TBST.

Detection: ECL Basic Kit (RM00020).

Exposure time: 90s.

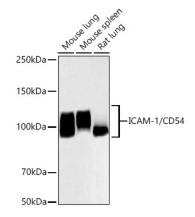

Western blot analysis of various lysates using ICAM-1/CD54 Rabbit PolymAb® (A26412-PM) at 1:5000 dilution incubated overnight at 4° C.


Secondary antibody: HRP-conjugated Goat anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.


Lysates/proteins: 25 μg per lane.


Blocking buffer: 3% nonfat dry milk in TBST.

Detection: ECL Basic Kit (RM00020).

Exposure time: 10s.

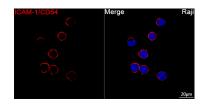

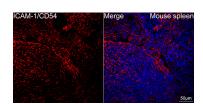

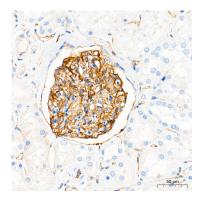

Confocal imaging of Raji cells using ICAM-1/CD54 Rabbit PolymAb® (A26412, dilution 1:4000) followed by a further incubation with Cy3 Goat Anti-Rabbit IgG (H+L) (AS007, dilution 1:500) (Red). DAPI was used for nuclear staining (Blue). Objective: 100x.

Confocal imaging of paraffin-embedded Mouse spleen tissue using ICAM-1/CD54 Rabbit PolymAb® (A26412, dilution 1:4000) followed by a further incubation with Cy3 Goat Anti-Rabbit IgG (H+L) (AS007, dilution 1:500) (Red). DAPI was used for nuclear staining (Blue). High pressure antigen retrieval performed with 0.01M Citrate Buffer (pH 6.0) prior to IF staining. Objective: 40x.


Immunohistochemistry analysis of paraffinembedded Human kidney tissue using ICAM-1/CD54 Rabbit PolymAb® (A26412-PM) at a dilution of 1:4000 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer (pH 9.0) prior to IHC staining.




Immunohistochemistry analysis of paraffinembedded Human lung tissue using ICAM-1/CD54 Rabbit PolymAb® (A26412-PM) at a dilution of 1:4000 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer (pH 9.0) prior to IHC staining.




Immunohistochemistry analysis of paraffinembedded Mouse liver tissue using ICAM-1/CD54 Rabbit PolymAb® (A26412-PM) at a dilution of 1:4000 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer (pH 9.0) prior to IHC staining.

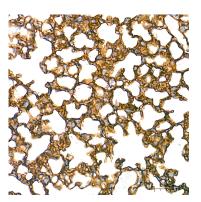

Immunohistochemistry analysis of paraffinembedded Mouse lung tissue using ICAM-1/CD54 Rabbit PolymAb® (A26412-PM) at a dilution of 1:4000 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer (pH 9.0) prior to IHC staining.

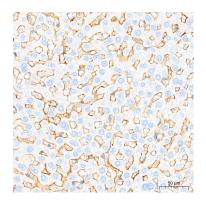

Immunohistochemistry analysis of paraffinembedded Rat liver tissue using ICAM-1/CD54 Rabbit PolymAb® (A26412-PM) at a dilution of 1:4000 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer (pH 9.0) prior to IHC staining.

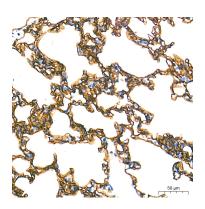

Immunohistochemistry analysis of paraffinembedded Rat lung tissue using ICAM-1/CD54 Rabbit PolymAb® (A26412-PM) at a dilution of 1:4000 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer (pH 9.0) prior to IHC staining