Ki67 Rabbit PolymAb®

Catalog No.: A26419



Basic Information

Observed MW

359kDa[[[

Calculated MW

351kDa

Category

Primary antibody

Applications

WB,IHC-P,IF/ICC,ELISA

Cross-Reactivity

Mouse, Rat

Background

Predicted to enable protein C-terminus binding activity. Involved in regulation of chromatin organization. Acts upstream of or within cholangiocyte proliferation; hepatocyte proliferation; and meiotic cell cycle. Located in chromosome; cytoplasm; and nucleolus. Is expressed in several structures, including alimentary system; brain; genitourinary system; sensory organ; and skin. Human ortholog(s) of this gene implicated in Crohn's disease; breast cancer; human immunodeficiency virus infectious disease; and pancreatic cancer. Orthologous to human MKI67 (marker of proliferation Ki-67).

Recommended Dilutions

WB 1:1000 - 1:5000

IHC-P 1:2000 - 1:12000

IF/ICC 1:2000 - 1:5000

ELISA Recommended starting

concentration is 1 µg/mL.

Please optimize the
concentration based on
your specific assay
requirements.

Immunogen Information

Gene IDSwiss Prot
17345
E9PVX6

Immunogen

Recombinant fusion protein containing a sequence corresponding to amino acids 1816-2163 of mouse Ki67 (NP_001074586.2).

Synonyms

Ki67; Ki-67; D630048A14Rik

Contact

a		400-999-6126
\bowtie		cn.market@abclonal.com.cn
\odot	Ī	www.abclonal.com.cn


Product Information

SourceIsotypePurificationRabbitIgGAffinity purification

Storage

Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS with 0.09% Sodium azide, 0.05% BSA, 50% glycerol, pH7.3.

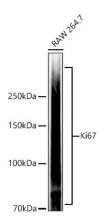

Western blot analysis of lysates from RAW 264.7 cells using Ki67 Rabbit PolymAb® (A26419) at 1:2000 dilution incubated overnight at 4° C.


Secondary antibody: HRP-conjugated Goat anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.


Lysates/proteins: 25 µg per lane.


Blocking buffer: 3% nonfat dry milk in TBST.

Detection: ECL Basic Kit (RM00020).


Exposure time: 5s.

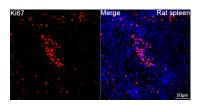

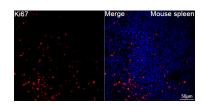


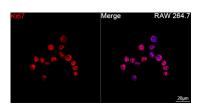


Confocal imaging of paraffin-embedded Rat spleen tissue using Ki67 Rabbit PolymAb® (A26419, dilution 1:4000) followed by a further incubation with Cy3 Goat Anti-Rabbit IgG (H+L) (AS007, dilution 1:500) (Red). DAPI was used for nuclear staining (Blue). High pressure antigen retrieval performed with 0.01M Citrate Buffer (pH 6.0) prior to IF staining. Objective: 40x.

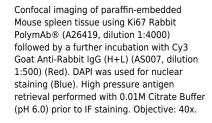


Confocal imaging of RAW 264.7 cells using Ki67 Rabbit PolymAb® (A26419, dilution 1:800) followed by a further incubation with Cy3 Goat Anti-Rabbit IgG (H+L) (AS007, dilution 1:500) (Red). DAPI was used for nuclear staining (Blue). Objective: 100x.

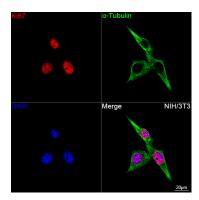


Confocal imaging of NIH/3T3 cells using Ki67 Rabbit PolymAb® (A26419, dilution 1:4000) followed by a further incubation with Cy3 Goat Anti-Rabbit IgG (H+L) (AS007, dilution 1:500) (Red). The cells were counterstained with α -Tubulin Mouse mAb (AC012, dilution 1:400) followed by incubation with ABflo® 488-conjugated Goat Anti-Mouse IgG (H+L) Ab (AS076, dilution 1:500) (Green). DAPI was used for nuclear staining (Blue). Objective: 100x.

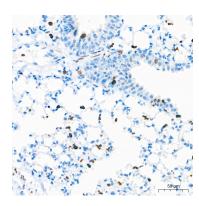


Immunohistochemistry analysis of paraffinembedded Mouse lung tissue using Ki67 Rabbit PolymAb® (A26419) at a dilution of 1:10000 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer(pH 9.0) prior to IHC staining.

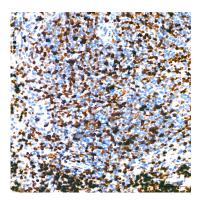

Immunohistochemistry analysis of paraffinembedded Mouse spleen tissue using Ki67 Rabbit PolymAb® (A26419) at a dilution of 1:10000 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer(pH 9.0) prior to IHC staining.

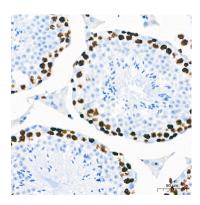

Immunohistochemistry analysis of paraffinembedded Mouse testis tissue using Ki67 Rabbit PolymAb® (A26419) at a dilution of 1:10000 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer(pH 9.0) prior to IHC staining.

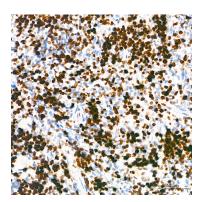

Immunohistochemistry analysis of paraffinembedded Rat spleen tissue using Ki67 Rabbit PolymAb® (A26419) at a dilution of 1:10000 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer(pH 9.0) prior to IHC staining.