


Validation Data

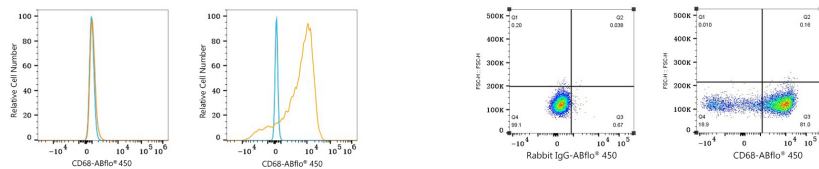

Flow cytometry: 1×10^6 NIH/3T3 cells (Low Expression, left) and RAW 264.7 cells (right) were intracellularly-stained with ABflo® 450 Rabbit anti-Mouse CD68 mAb (A26451, 5 μ l/Test, orange line) or ABflo® 450 Rabbit IgG isotype control (5 μ l/Test, blue line). Non-fluorescently stained cells were used as blank control (red line).

Flow cytometry: 1×10^6 RAW 264.7 cells were intracellularly-stained with ABflo® 450 Rabbit IgG isotype control (5 μ l/Test, left) or ABflo® 450 Rabbit anti-Mouse CD68 mAb (A26451, 5 μ l/Test, right).