Leader in Biomolecular Solutions for Life Science

Calpain 2 Rabbit pAb

Catalog No.: A2671



Basic Information

Observed MW 80kDa

Calculated MW 80kDa

Category Primary antibody

Applications WB,IHC-P,ELISA

Cross-Reactivity Human, Mouse, Rat

Background

The calpains, calcium-activated neutral proteases, are nonlysosomal, intracellular cysteine proteases. The mammalian calpains include ubiquitous, stomach-specific, and muscle-specific proteins. The ubiquitous enzymes consist of heterodimers with distinct large, catalytic subunits associated with a common small, regulatory subunit. This gene encodes the large subunit of the ubiquitous enzyme, calpain 2. Multiple heterogeneous transcriptional start sites in the 5' UTR have been reported. Two transcript variants encoding different isoforms have been found for this gene.

Recommended Dilutions

WB	1:500 - 1:1000
IHC-P	1:50 - 1:200
ELISA	Recommended starting concentration is 1 µg/mL. Please optimize the concentration based on your specific assay requirements.

Immunogen Information

Gene ID 824

Swiss Prot P17655

Immunogen

Recombinant protein (or fragment). This information is considered to be commercially sensitive.

Synonyms

CANP2; mCANP; CANPL2; CANPml; Calpain 2

Contact


6	400-999-6126
\mathbf{X}	cn.market@abclonal.com.cn
€	www.abclonal.com.cn
€	www.abclonal.com.c


Product Information

Source Rabbit **Isotype** IgG **Purification** Affinity purification


Storage

Store at -20°C. Avoid freeze / thaw cycles. Buffer: PBS with 0.09% Sodium azide,50% glycerol,pH7.3.

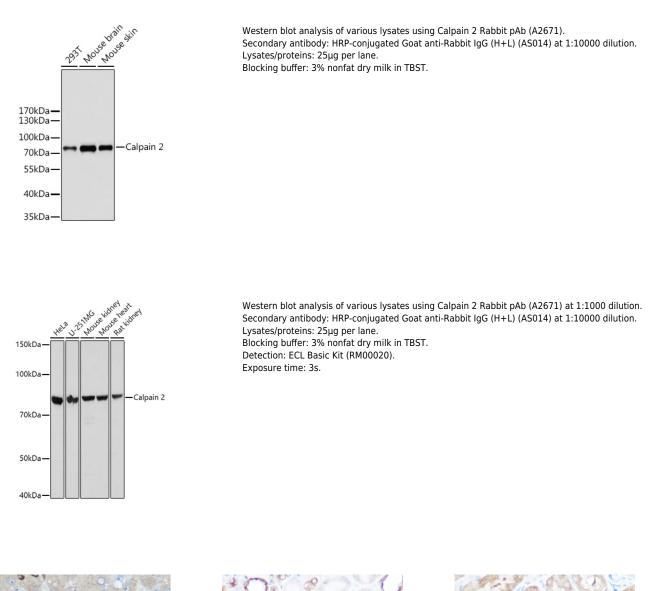


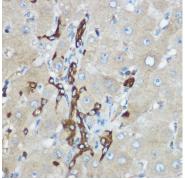

Immunohistochemistry analysis of paraffinembedded Human liver using Calpain 2 Rabbit pAb (A2671) at dilution of 1:100 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate buffer (pH 6.0) prior to IHC staining.

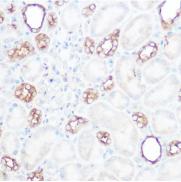

Immunohistochemistry analysis of paraffinembedded Mouse kidney using Calpain 2 Rabbit pAb (A2671) at dilution of 1:100 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate buffer (pH 6.0) prior to IHC staining.

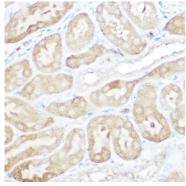

Immunohistochemistry analysis of paraffinembedded Rat kidney using Calpain 2 Rabbit pAb (A2671) at dilution of 1:100 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate buffer (pH 6.0) prior to IHC staining.