[KO Validated] Ki67 Rabbit PolymAb®

Catalog No.: A26755PM KO Validated



Basic Information

Observed MW 359kDa/359kDa

Calculated MW 350kDa/359kDa

Category Primary antibody

Applications WB,IHC-P,IF/ICC,ELISA

Cross-Reactivity Human, Mouse, Rat

Background

Enables protein C-terminus binding activity. Involved in regulation of chromosome segregation and regulation of mitotic nuclear division. Located in chromosome; nuclear body; and nucleolus. Colocalizes with condensed chromosome. Implicated in Crohn's disease; breast cancer; human immunodeficiency virus infectious disease; and pancreatic cancer. Biomarker of several diseases, including Barrett's esophagus; autoimmune disease of musculoskeletal system (multiple); endocrine gland cancer (multiple); gastrointestinal system cancer (multiple); and interstitial cystitis.

Recommended Dilutions

WB	1:1000 - 1:6000
IHC-P	1:1000 - 1:8000
IF/ICC	1:200 - 1:800
ELISA	Recommended starting concentration is 1 µg/mL. Please optimize the concentration based on your specific assay requirements.

Immunogen Information

Gene ID 17345 Swiss Prot E9PVX6

Immunogen

Synthetic peptide. This information is considered to be commercially sensitive.

Synonyms KIA; MIB-; MIB-1; PPP1R105

Contact


6	400-999-6126
\bowtie	cn.market@abclonal.com.cn
€	www.abclonal.com.cn


Product Information

Source Rabbit **lsotype** IgG Purification Affinity purification


Storage

Store at -20°C. Avoid freeze / thaw cycles. Buffer: PBS with 0.09% Sodium azide,0.05% BSA,50% glycerol,pH7.3.

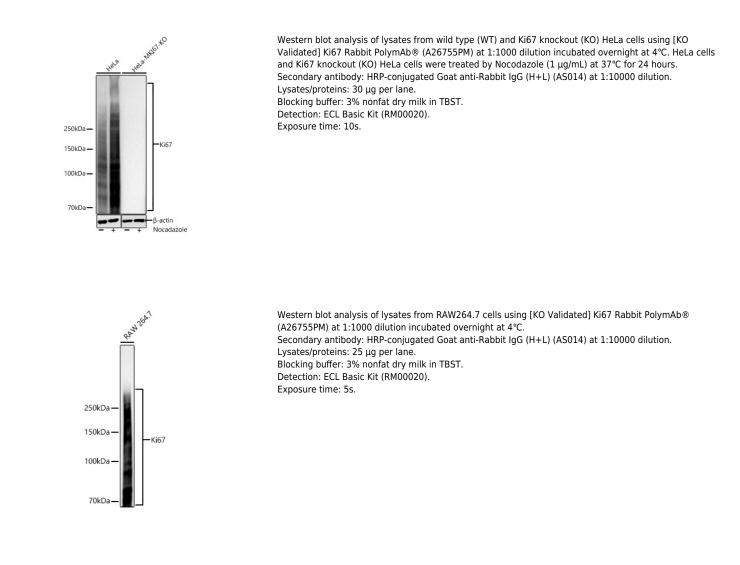


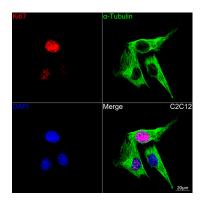

Confocal imaging of C2C12 cells using [KO Validated] Ki67 Rabbit PolymAb® (A26755PM, dilution 1:200) followed by a further incubation with Cy3 Goat Anti-Rabbit IgG (H+L) (AS007, dilution 1:500) (Red). The cells were counterstained with α -Tubulin Mouse mAb (AC012, dilution 1:400) followed by incubation with ABflo® 488-conjugated Goat Anti-Mouse IgG (H+L) Ab (AS076, dilution 1:500) (Green). DAPI was used for nuclear staining (Blue). Objective: 100x.

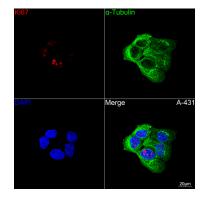


Confocal imaging of A-431 cells using [KO Validated] Ki67 Rabbit PolymAb® (A26755PM, dilution 1:200) followed by a further incubation with Cy3 Goat Anti-Rabbit IgG (H+L) (AS007, dilution 1:500) (Red). The cells were counterstained with α -Tubulin Mouse mAb (AC012, dilution 1:400) followed by incubation with ABflo® 488-conjugated Goat Anti-Mouse IgG (H+L) Ab (AS076, dilution 1:500) (Green). DAPI was used for nuclear staining (Blue). Objective: 100x.

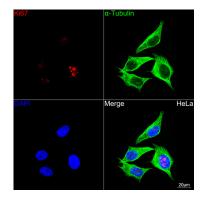


Confocal imaging of HeLa cells using [KO Validated] Ki67 Rabbit PolymAb® (A26755PM, dilution 1:200) followed by a further incubation with Cy3 Goat Anti-Rabbit IgG (H+L) (AS007, dilution 1:500) (Red). The cells were counterstained with α -Tubulin Mouse mAb (AC012, dilution 1:400) followed by incubation with ABflo® 488-conjugated Goat Anti-Mouse IgG (H+L) Ab (AS076, dilution 1:500) (Green). DAPI was used for nuclear staining (Blue). Objective: 100x.


Validation Data

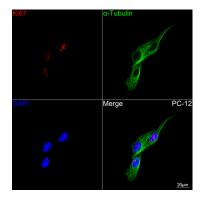


Confocal imaging of PC-12 cells using [KO Validated] Ki67 Rabbit PolymAb® (A26755PM, dilution 1:200) followed by a further incubation with Cy3 Goat Anti-Rabbit IgG (H+L) (AS007, dilution 1:500) (Red). The cells were counterstained with α -Tubulin Mouse mAb (AC012, dilution 1:400) followed by incubation with ABflo® 488-conjugated Goat Anti-Mouse IgG (H+L) Ab (AS076, dilution 1:500) (Green). DAPI was used for nuclear staining (Blue). Objective: 100x.

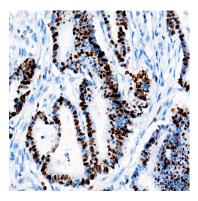


Immunohistochemistry analysis of paraffinembedded Human colon carcinoma tissue using [KO Validated] Ki67 Rabbit PolymAb® (A26755PM) at a dilution of 1:1000 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer (pH 9.0) prior to IHC staining.

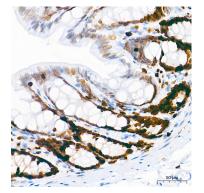


Immunohistochemistry analysis of paraffinembedded Rat colon tissue using [KO Validated] Ki67 Rabbit PolymAb® (A26755PM) at a dilution of 1:1000 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer (pH

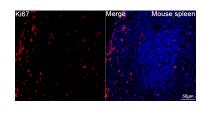


Confocal imaging of paraffin-embedded Mouse spleen tissue using [KO Validated] Ki67 Rabbit PolymAb® (A26755PM, dilution 1:200) followed by a further incubation with Cy3 Goat Anti-Rabbit IgG (H+L) (AS007, dilution 1:500) (Red). DAPI was used for nuclear staining (Blue). High pressure antigen retrieval performed with 0.01M Citrate Buffer (pH 6.0) prior to IF staining. Objective: 40x.

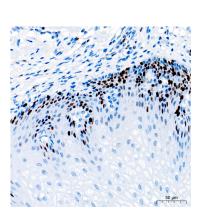


Immunohistochemistry analysis of paraffinembedded Human esophagus tissue using [KO Validated] Ki67 Rabbit PolymAb® (A26755PM) at a dilution of 1:1000 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer (pH 9.0) prior to IHC staining.

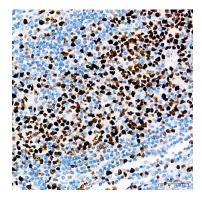

Immunohistochemistry analysis of paraffinembedded Human tonsil tissue using [KO Validated] Ki67 Rabbit PolymAb® (A26755PM) at a dilution of 1:1000 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer (pH 9.0) prior to IHC staining.

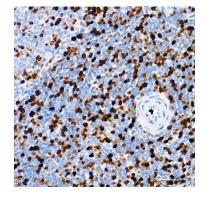

Immunohistochemistry analysis of paraffinembedded Rat spleen tissue using [KO Validated] Ki67 Rabbit PolymAb® (A26755PM) at a dilution of 1:1000 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer (pH

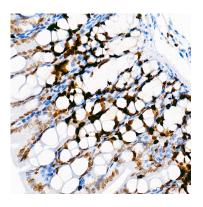

Immunohistochemistry analysis of paraffinembedded Mouse colon tissue using [KO Validated] Ki67 Rabbit PolymAb® (A26755PM) at a dilution of 1:1000 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer (pH 9.0) prior to IHC staining.

9.0) prior to IHC staining.

9.0) prior to IHC staining.