FOSL2 Rabbit pAb

Catalog No.: A2729 1 Publications



Basic Information

Observed MW

35kDa

Calculated MW

35kDa

Category

Primary antibody

Applications

WB,IHC-P,IF/ICC,ELISA

Cross-Reactivity

Human, Rat

Background

The Fos gene family consists of 4 members: FOS, FOSB, FOSL1, and FOSL2. These genes encode leucine zipper proteins that can dimerize with proteins of the JUN family, thereby forming the transcription factor complex AP-1. As such, the FOS proteins have been implicated as regulators of cell proliferation, differentiation, and transformation.

Recommended Dilutions

WB 1:500 - 1:2000

IHC-P 1:100 - 1:200

IF/ICC 1:50 - 1:200

ELISA Recommended starting

concentration is 1 µg/mL.

Please optimize the
concentration based on
your specific assay
requirements.

Immunogen Information

Gene ID2355

Swiss Prot
P15408

Immunogen

Recombinant fusion protein containing a sequence corresponding to amino acids 187-326 of human FOSL2 (NP_005244.1).

Synonyms

FRA2; FOSL2

Contact

	400-999-6126
×	cn.market@abclonal.com.cn
	www.abclonal.com.cn


Product Information

SourceIsotypePurificationRabbitIgGAffinity purification

Storage

Store at -20 $^{\circ}\text{C}.$ Avoid freeze / thaw cycles.

Buffer: PBS with 0.02% sodium azide,50% glycerol,pH7.3.

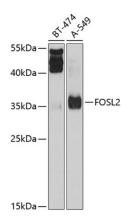


Western blot analysis of various lysates using FOSL2 Rabbit pAb (A2729) at 1:1000 dilution. Secondary antibody: HRP-conjugated Goat anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.


Lysates/proteins: 25µg per lane.

Blocking buffer: 3% nonfat dry milk in TBST. Detection: ECL Basic Kit (RM00020).

Exposure time: 30s.

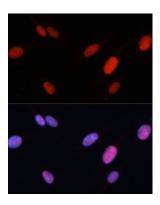

Immunofluorescence analysis of U-2 OS cells using FOSL2 Rabbit pAb (A2729) at dilution of 1:100. Secondary antibody: Cy3conjugated Goat anti-Rabbit IgG (H+L) (AS007) at 1:500 dilution. Blue: DAPI for nuclear staining.

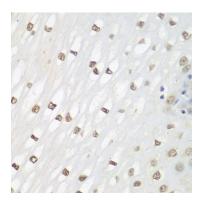

Immunohistochemistry analysis of paraffinembedded Human esophagus using FOSL2 Rabbit pAb (A2729) at dilution of 1:100 (40x lens). Microwave antigen retrieval performed with 0.01M PBS Buffer (pH 7.2) prior to IHC staining.