SUMO4 Rabbit pAb

Catalog No.: A3100



Basic Information

Observed MW

18kDa

Calculated MW

11kDa

Category

Primary antibody

Applications

ELISA,WB,IHC-P

Cross-Reactivity

Human, Mouse

Background

This gene is a member of the SUMO gene family. This family of genes encode small ubiquitin-related modifiers that are attached to proteins and control the target proteins' subcellular localization, stability, or activity. The protein described in this record is located in the cytoplasm and specifically modifies IKBA, leading to negative regulation of NF-kappa-B-dependent transcription of the IL12B gene. A specific polymorphism in this SUMO gene, which leads to the M55V substitution, has been associated with type I diabetes. The RefSeq contains this polymorphism.

Recommended Dilutions

WB 1:1000 - 1:2000

IHC-P 1:100 - 1:300

Immunogen Information

Gene ID387082

Swiss Prot
Q6EEV6

Immunogen

A synthetic peptide corresponding to a sequence within amino acids 40-95 of human SUMO4 (Q6EEV6).

Synonyms

IDDM5; SMT3H4; SUMO-4; dJ281H8.4; SUMO4

Contact

a	400-999-6126
×	cn.market@abclonal.com.cn
\overline{a}	www.abclonal.com.cn


Product Information

SourceIsotypePurificationRabbitIgGAffinity purification


Storage

Store at -20°C. Avoid freeze / thaw cycles.


Buffer: PBS with 0.02% sodium azide,50% glycerol,pH7.3.

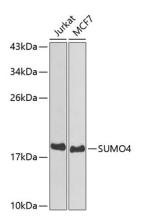


Western blot analysis of various lysates using SUMO4 Rabbit pAb (A3100) at 1:800 dilution. Secondary antibody: HRP Goat Anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution. Lysates/proteins: $25\mu g$ per lane. Blocking buffer: 3% nonfat dry milk in TBST.

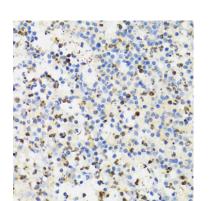

Immunohistochemistry analysis of paraffinembedded human appendicitis using SUMO4 Rabbit pAb (A3100) at dilution of 1:100 (40x lens).Perform microwave antigen retrieval with 10 mM PBS buffer pH 7.2 before commencing with IHC staining protocol.

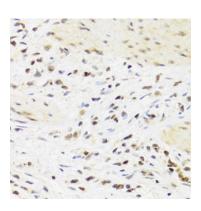

Immunohistochemistry analysis of paraffinembedded human breast cancer using SUMO4 Rabbit pAb (A3100) at dilution of 1:100 (40x lens).Perform microwave antigen retrieval with 10 mM PBS buffer pH 7.2 before commencing with IHC staining protocol.

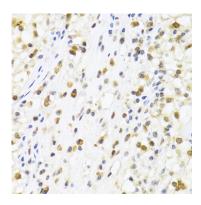

Immunohistochemistry analysis of paraffinembedded human kidney cancer using SUMO4 Rabbit pAb (A3100) at dilution of 1:100 (40x lens).Perform microwave antigen retrieval with 10 mM PBS buffer pH 7.2 before commencing with IHC staining protocol.