CD35/CR1 Rabbit mAb

Catalog No.: A3661 Recombinant



Basic Information

Observed MW

Calculated MW

224kDa

Category

Primary antibody

Applications

ELISA,IHC-P,IF/ICC

Cross-Reactivity

Human, Mouse, Rat

CloneNo number

ARC2065

Background

This gene is a member of the receptors of complement activation (RCA) family and is located in the 'cluster RCA' region of chromosome 1. The genome is polymorphic at this locus with allele-specific splice variants encoding different isoforms, based on the presence/absence of long homologous repeats (LHRs). The gene encodes a monomeric single-pass type I membrane glycoprotein found on erythrocytes, leukocytes, glomerular podocytes, and splenic follicular dendritic cells. The Knops blood group system is a system of antigens located on this protein. The protein mediates cellular binding to particles and immune complexes that have activated complement. Decreases in expression of this protein and/or mutations in this gene have been associated with gallbladder carcinomas, mesangiocapillary glomerulonephritis, systemic lupus erythematosus, sarcoidosis and Alzheimer's disease. Mutations in this gene have also been associated with a reduction in Plasmodium falciparum rosetting, conferring protection against severe malaria.

Recommended Dilutions

IHC-P 1:50 - 1:200

IF/ICC 1:50 - 1:200

Immunogen Information

Gene ID Swiss Prot 1378 P17927

Immunogen

A synthetic peptide corresponding to a sequence within amino acids 1940-2039 of human CD35/CR1 (P17927).

Synonyms

KN; C3BR; C4BR; CD35; CD35/CR1

Contact

6	400-999-6126
\bowtie	cn.market@abclonal.com.cn
•	www.abclonal.com.cn


Product Information

SourceIsotypePurificationRabbitIgGAffinity purification


Storage

Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS with 0.02% sodium azide, 0.05% BSA, 50% glycerol, pH7.3.

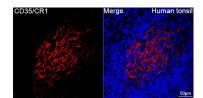

Confocal imaging of paraffin-embedded Human tonsil tissue using CD35/CR1 Rabbit mAb (A3661, dilution 1:100) followed by a further incubation with Cy3 Goat Anti-Rabbit IgG (H+L) (AS007, dilution 1:500) (Red). DAPI was used for nuclear staining (Blue). Objective: 40x. Perform high pressure antigen retrieval with 0.01 M citrate buffer (pH 6.0) prior to IF staining.

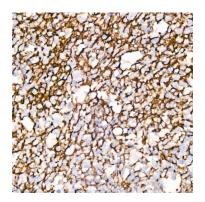

Immunohistochemistry analysis of CD35/CR1 in paraffin-embedded human tonsil using CD35/CR1 Rabbit mAb (A3661) at dilution of 1:50 (40x lens).Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.