Coagulation Protein C Rabbit mAb

Catalog No.: A4532 Recombinant



Basic Information

Observed MW

52kDa

Calculated MW

52kDa

Category

Primary antibody

Applications

ELISA,WB,IF/ICC

Cross-Reactivity

Human, Mouse, Rat

CloneNo number

ARC1023

Background

This gene encodes a vitamin K-dependent plasma glycoprotein. The encoded protein is cleaved to its activated form by the thrombin-thrombomodulin complex. This activated form contains a serine protease domain and functions in degradation of the activated forms of coagulation factors V and VIII. Mutations in this gene have been associated with thrombophilia due to protein C deficiency, neonatal purpura fulminans, and recurrent venous thrombosis.

Recommended Dilutions

WB 1:500 - 1:2000

IF/ICC 1:50 - 1:200

Immunogen Information

Gene ID Swiss Prot 5624 P04070

Immunogen

Recombinant fusion protein containing a sequence corresponding to amino acids 300-461 of human Coagulation Coagulation Protein C (P04070).

Synonyms

PC; APC; PROC1; THPH3; THPH4; Coagulation Protein C

Contact

| a | 400-999-6126 |
|----------|---------------------------|
| \sim | cn.market@abclonal.com.cn |
| ⊙ | www.abclonal.com.cn |

Product Information


SourceIsotypePurificationRabbitIgGAffinity purification

Storage

Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS with 0.02% sodium azide, 0.05% BSA, 50% glycerol, pH7.3.

Validation Data

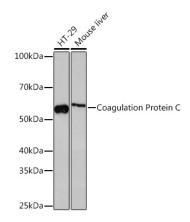

Western blot analysis of various lysates using Coagulation Coagulation Protein C Rabbit mAb (A4532) at 1:1000 dilution.


Secondary antibody: HRP Goat Anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.

Lysates/proteins: 25µg per lane.

Blocking buffer: 3% nonfat dry milk in TBST.

Detection: ECL Basic Kit (RM00020).

Exposure time: 3min.

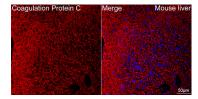

Confocal imaging of paraffin-embedded Mouse liver tissue using Coagulation Protein C Rabbit mAb (A4532, dilution 1:100) followed by a further incubation with Cy3 Goat Anti-Rabbit IgG (H+L) (AS007, dilution 1:500) (Red). DAPI was used for nuclear staining (Blue). Objective: 40x. Perform high pressure antigen retrieval with 0.01 M citrate buffer (pH 6.0) prior to IF staining.